Betriebsanleitung SensorReady 3D

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The 30W interchangeable AC clip adaptor series is the latest offering from a world leader in charging adapter technology and manufacturing - Pihong Technology. Pihong 30W interchangeable AC clip adaptor series are designed to work with most 12V/15V/24V/56V powered devices and conform to most of the industry standards and meet DoE level VI energy efficiency requirements. We ensure our products meet or exceed all the industry performance and safety requirements.

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
30W Interchangeable AC Clip Adapter

PSAA30R-XXX-R Series

● 12V/15V/24V/56V ● DoE 6 ● Interchangeable AC Clip ● Cable ● Black

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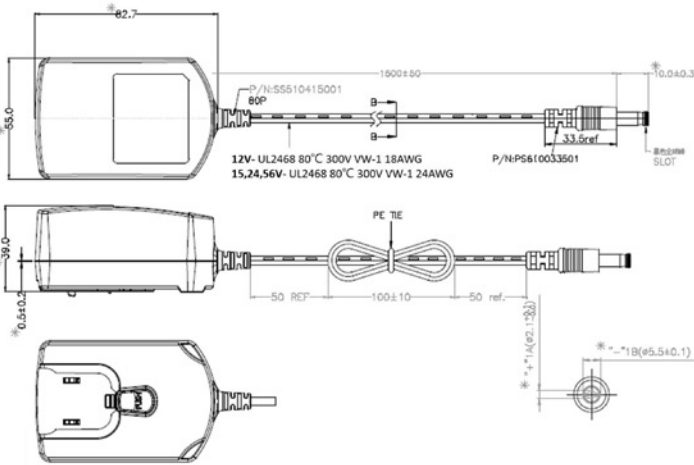declares that the following product

Product name: **SensorReady 3D**
 Article no.: **00-76-901-4830**
 Product type: **USB interface for incremental
 length measuring systems**

is in conformity with the following European Union harmonisation legislation, when operated in accordance with the operating manual:

EMC directive 2014/30/EU
Low voltage directive 2014/35/EU
RoHS directive 2011/65/EU

The following harmonised standards have been applied:

DIN EN 61010-1:2011-07
DIN EN 61326-1:2013-07

In case of a modification of the product, which was not approved by the manufacturer, this Declaration of Conformity loses its validity.

This Declaration of Conformity is issued under the sole responsibility of the manufacturer.

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Approved by  Uwe Schnitter | 13.09.2016

Issued by Christopher Koch | 13.09.2016 | 10:30

Proved to be QMS conform by  Christopher Koch | 13.09.2016

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DE

Table of Contents

1	Basics of the Operating Manual	61
1.1	Symbols and Abbreviations	62
2	Delivery Contents	63
2.1	Optional Accessories	63
3	Product Description	64
3.1	Identification	64
3.2	Area of Application	64
3.3	Compliance with Directives and Standards	65
3.4	Structure and Ports	66
4	Safety	67
4.1	Intended Use	67
4.2	Installation Requirements	67
4.3	General Safety Notes	68
4.4	Warning Note Categories	70
4.5	Hazard Symbols Used	70
5	Putting into Service	71
5.1	Putting SensorReady 3D into Service	71
5.2	Establishing Connection to a PC	75
6	Operation	77
6.1	Operation with GetPos	77
6.2	Operation via Instruction Set	85
6.3	Using Position Correction	86
7	Cleaning	88
8	Maintenance	88
9	Troubleshooting	89
10	Service	92
10.1	Replacement Components	93
11	Taking out of Service	94
12	Disposal	95

- 13 Technical Data 96
 - 13.1 General Information 96
 - 13.2 Encoder Interface 98
 - 13.3 Ports 99
 - 13.4 External 24 V/30 W Power Supply (PSAA30R-240) 106
- 14 Declaration of Conformity 108
- 15 Index 109

1 Basics of the Operating Manual

This operating manual gives you the information which is needed for trouble-free and safe use of the SensorReady 3D. You must have fully read and understood this operating manual and especially the safety instructions before you put into service, operate, clean or dispose of the SensorReady 3D.

Märzhäuser Sensotech accepts no liability for damage and operating faults which result from ignoring this operating manual.

This operating manual is applicable for a SensorReady 3D with firmware version 1.02 and higher.

1.1 Symbols and Abbreviations

Symbol	Meaning
	Additional information and useful notes
	Start of an instruction
1, 2, 3 ...	Steps within one instruction
»	Cross-reference to another part of this operating manual, another document or a website

Abbreviation	Meaning
Fig.	Figure
approx.	approximately, roughly
if nec.	if necessary
max.	maximum
min.	minimum
typ.	typically
e.g.	for example

2 Delivery Contents

No.	Designation	Quantity
1	SensorReady 3D	1
2	MR encoder cable	1
3	USB interface cable	1
4	External power supply 24 V/30 W incl. 4 socket adapters (EU/GB/USA/AUS)	1
5	Operating manual	1
6	Installation CD incl. driver, operating manual, software <i>GetPos</i>	1

2.1 Optional Accessories

The following optional accessories are available for the SensorReady 3D:

Designation	Order no.
Customised encoder cables	On request
Mounting plate for top hat rail installation	On request

For questions about ordering, please contact our sales team:

Tel.: **+49 6441 9116-10**
 Fax: **+49 6441 9116-40**
 E-mail: **sales@marzhauser.com**

If you want to make your own cables, you can find the pin assignments of all the SensorReady 3D ports here:

» *Chapter 13.3 Ports, page 99*

3 Product Description

3.1 Identification

You can find all required information for clear identification of the SensorReady 3D on the nameplate. The nameplate is attached to the back of the SensorReady 3D.

3.2 Area of Application

The SensorReady 3D is used for signal evaluation during the measurement of objects by means of positioning systems such as manual measuring stages or Z-axes which are equipped with measuring systems.

Using the integrated position correction, systematic length errors of the positioning system can be corrected. Evaluating reference signals creates an absolute reference of the position correction values to the machine zero point.

Operation of the SensorReady 3D takes place via a USB interface on a PC. Evaluation via the PC is carried out with the *GetPos* software or directly by using the instruction set. The SensorReady 3D can also be programmed using the instruction set.

The SensorReady 3D is compatible with all positioning systems from Märzhäuser Wetzlar that have integrated measuring systems.

Here you can find further information on the magnetic measuring system from Märzhäuser Sensotech:

» www.marzhauser-st.com

3.3 Compliance with Directives and Standards


Märzhäuser Sensotech hereby declares that the SensorReady 3D was designed and produced in compliance with the EU directives and harmonised standards listed in the Declaration of Conformity.

The Declaration of Conformity for the SensorReady 3D can be found here:

» *Chapter 14 Declaration of Conformity, page 108*

3.4 Structure and Ports

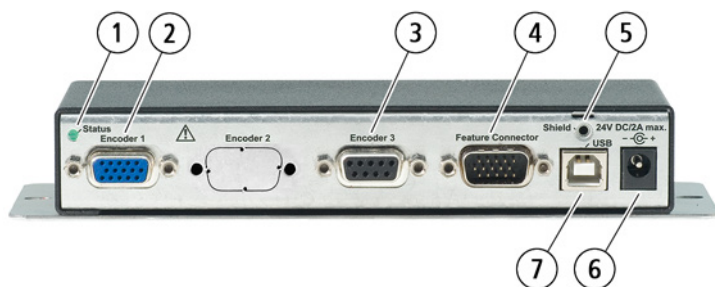No.	Designation	Function
1	Status	Status indicator
2	Encoder 1	Port for measuring systems 1 and 2 (<i>MR/analog 5Vss</i>)
3	Encoder 3	Port for measuring system 3 (<i>MR/ analog 5Vss, 1Vss, TTL (RS-422)</i>)
4	Feature Connector	Optional functions (e.g. reference switch input 3rd measuring system, RS-232 interface)
5	Shield	Shield port/ground port
6	24V DC/ 2A max.	Power supply port (24 V)
7	USB	USB port (PC interface)

4 Safety

4.1 Intended Use

The SensorReady 3D is intended and designed for the connection and evaluation of

- 2 incremental measuring systems of the type *MR/analog 5Vss* and
- 1 further incremental measuring system without external power supply of the type *MR/analog 5Vss, 1Vss* or *TTL (RS-422)*.

On request, your SensorReady 3D can also be upgraded for the connection of an incremental measuring system of the type *TTL (non-differential)* to port *Encoder 3*. Please contact our Customer Service Center if necessary:

» *Chapter 10 Service, page 92*

4.2 Installation Requirements

Category	Condition
Installation location	Closed rooms, no risk of explosion
Installation height	Max. 2,000 m
Environmental temperature	+5...+40 °C
Storage temperature	-20...+70 °C
Cooling	Normal convection
Humidity	Max. 85% non-condensing

4.3 General Safety Notes

- Only use the SensorReady 3D for the purposes which are described in this operating manual.
- Read and make sure you understand the operating manual for the measuring systems and peripherals you would like to connect the SensorReady 3D to.
- Only connect measuring systems from Märzhäuser Wetzlar to port *Encoder 1*.
- Only connect incremental measuring systems of the type *MR/analog 5Vss, 1Vss or TTL (RS-422)* without external power supply to port *Encoder 3*.
- Only use shielded cables for connecting measuring systems, PCs and other peripherals.
- The SensorReady 3D is insulated to the power supply. Ensure that no dangerous voltages or currents are induced by peripheral devices in the SensorReady 3D.
- Only use the 24 V/30 W power supply provided to power the SensorReady 3D.
- Ensure that your mains voltage matches the label on the power supply before connecting the SensorReady 3D to the mains.
- Disconnect the SensorReady 3D from the mains before connecting measuring systems, a PC or other peripherals to the SensorReady 3D or disconnecting them from the SensorReady 3D.
- Set the SensorReady 3D up in such a way that you can immediately disconnect the power supply from the mains, without hindrance, in case of emergency.
- Ensure that no chips, liquids, objects or substances enter the interior of the SensorReady 3D or the power supply. Protect the SensorReady 3D from dust and dirt as far as possible.

- The SensorReady 3D is not protected against corrosive, infectious, radioactive or other substances which are harmful to health. When handling such substances, ensure that all legal requirements are fulfilled, especially the national regulations on accident prevention.
- Immediately replace the power supply if the power supply housing or cable show signs of damage.
- Disconnect the SensorReady 3D from the mains immediately if there are malfunctions, defects or safety flaws. Secure the SensorReady 3D from being reconnected to the mains.
- Never open the housing of the SensorReady 3D or the power supply provided. In case of violation, all warranty and liability claims will become void.
- Only have repair work carried out by our Customer Service Center.
 - » *Chapter 10 Service, page 92*
- Ensure that the SensorReady 3D, in conjunction with your application, meets the applicable safety regulations and legal requirements.
- Store this operating manual so that it is readily available at the installation location.
- If you lend or sell the SensorReady 3D you must include all the documents provided.

4.4 Warning Note Categories

Warning notes are divided into the following categories:

WARNING

Hazardous situation that can cause serious injury or death if not avoided.

NOTICE

Warning of possible damage or malfunction

4.5 Hazard Symbols Used

The following hazard symbols can be used on the SensorReady 3D and/or in warning notes in this operating manual:

Symbol	Meaning
	Warning of a general hazard which can be avoided by observing the operating manual
	Warning of hazards caused by electric current
	Warning of hazards caused by substances which are harmful to health

5 Putting into Service

5.1 Putting SensorReady 3D into Service

5.1.1 Connecting Measuring System 1 and 2

NOTICE

Possible equipment damage through putting into/taking out of service with existing mains power connection.

Disconnect the SensorReady 3D from the mains before connecting/disconnecting measuring systems, a PC or other peripherals to/from the SensorReady 3D.

Incorrect data transmission through unshielded cables possible.

Only use shielded cables for connecting measuring systems, PCs and other peripherals.

► **Putting SensorReady 3D into service**

- 1 Read and make sure you understand the *Safety* chapter before you start set-up and putting into service.
» *Chapter 4 Safety, page 67*
- 2 Read and make sure you understand the operating manual of the peripherals you would like to connect the SensorReady 3D to.
- 3 Ensure that the SensorReady 3D has a secure footing and is mounted stably.
- 4 Ensure that the SensorReady 3D is not connected to the mains.

NOTICE

Equipment damage possible by connecting inappropriate measuring systems.

Only connect measuring systems from Märzhäuser Wetzlar to port *Encoder 1*.

- 5 Connect your measuring systems 1 and 2 to port *Encoder 1* on the SensorReady 3D. Use the MR encoder cable for this.

Note: The reference switches for your measuring systems 1 and 2 are thus also connected to the SensorReady 3D.

WARNING

Possible electric shock and equipment damage by connection to external or damaged power supplies.

Only use the 24 V/30 W power supply provided. Immediately replace the power supply if the power supply housing or cable show signs of damage. Only use replacement components from Märzhäuser Sensotech.

Possible electric shock and equipment damage by connection to a non-specified mains.

Ensure that the voltage of your mains matches the label on the power supply before connecting the power supply to the mains.

- 6 Connect the power supply cable to the port *24V DC/2A max.*
- 7 Connect the power supply to the mains.

The status LED lights up.

- 8 Move your positioning system at least 1 cm in both directions.

Note: The measuring system is automatically calibrated in this way. The calibration is kept after switching off the SensorReady 3D and is continually updated during operation.

Putting into service is now complete. You can now perform measurements with the SensorReady 3D.

5.1.2 Connecting Measuring System 3

You can use an additional 3rd measuring system (probe, linear encoder or rotary encoder) with the SensorReady 3D to supplement measuring systems 1 and 2.

At the port *Encoder 3*, you can connect an incremental measuring system without external power supply from Märzhäuser Sensotech (*MR/analog 5Vss*) or another manufacturer to an interface of type *1Vss* and *TTL (RS-422)*.

You can get suitable encoder cables for your 3rd measuring system as an optional accessory from our Customer Service Center:

» *Chapter 10 Service, page 92*

► Connecting 3rd measuring system

- 1 Read and make sure you understand the operating manual of your 3rd measuring system.

NOTICE

Possible equipment damage through putting into/taking out of service with existing mains power connection.

Disconnect the SensorReady 3D from the mains before connecting/disconnecting measuring systems, a PC or other peripherals to/from the SensorReady 3D.

- 2 Disconnect the power supply of the SensorReady 3D from the mains.

NOTICE

Equipment damage possible from by connecting inappropriate measuring systems.

Only connect incremental measuring systems of type *MR/analog 5Vss*, *1Vss* or *TTL (RS-422)* without external power supply to port *Encoder 3*.

Incorrect data transmission through unshielded cables possible.

Only use shielded cables for connecting measuring systems, PCs and other peripherals.

3 Connect your 3rd measuring system to port *Encoder 3* on the SensorReady 3D. Use an appropriate encoder cable for this.

4 If you are using a reference switch in connection with your 3rd measuring system, connect the reference switch to the port *Feature Connector*. Use an appropriate cable for this.

» *Chapter 2.1 Optional Accessories, page 63*

Note: If your 3rd measuring system has reference marks, the reference signal is transmitted via port *Encoder 3* to the SensorReady 3D.

5 Connect the power supply of the SensorReady 3D to the mains.

The status LED lights up. The SensorReady 3D is ready for operation.

6 If you are using a measuring system of type *MR/analog 5Vss*, you must calibrate your 3rd measuring system. For this purpose, move your third axis at least 1 cm in both directions.

Note: The calibration is kept after switching off the SensorReady 3D and is continually updated during operation. Measuring systems of type *1Vss* and *TTL (RS-422)* do not need to be calibrated.

Your 3rd measuring system can now be evaluated by the SensorReady 3D.

5.2 Establishing Connection to a PC

Communication between SensorReady 3D and PC takes place via a virtual COM interface. To connect the SensorReady 3D to your PC via USB cable, you require *Windows 2000*, *Windows XP (32/64 Bit)*, *Windows Vista (32/64 Bit)*, *Windows 7 (32/64 Bit)*, *Windows 8 (32/64 Bit)* or *Windows 10 (32/64 Bit)*.

Observe the following notes, regardless of the version of your operating system:

- Allow the drivers to install automatically from the operating system.
- Do not connect to *Windows Update*.
- Continue installing the drivers despite the message about the missing *Windows-Logo-Test*.

► Establishing connection to a PC using a USB cable

- 1 Disconnect the power supply of the SensorReady 3D from the mains.
- 2 Switch the PC on.
- 3 Insert the installation CD in the CD drive.
- 4 Connect the *USB* port on the SensorReady 3D to a USB port on your PC. Use the USB interface cable for this.
- 5 Connect the power supply of the SensorReady 3D to the mains.

The operating system automatically recognises a new device (SensorReady 3D).

- 6 Follow the installation instructions of your operating system.

Note: If the driver for the SensorReady 3D is not automatically found on the installation CD, enter the driver storage location manually (INF file). You can find the drivers for the SensorReady 3D on the installation CD in the `..\Drivers` folder.

Once installation is completed, the SensorReady 3D is connected to your PC.

► Checking the version of the installed driver

If necessary, you can check the version of the installed driver.

- 1 Open the device manager:
 - **Windows 2000, Windows XP:** *Start → Control panel → System → Hardware → Device manager*
 - **Windows Vista:** *Start → right-click on Computer → Device manager*
 - **Windows 7:** *Start → Control panel → System → Device manager*
 - **Windows 8:** Right-click in the left, lower corner of the desktops/start screen → *Device manager*
 - **Windows 10:** Right-click the Windows symbol → *Device manager*
- 2 Expand the column *Ports (COM and LPT)* by double-clicking.
- 3 Right-click the item *SensorReady 3D* to open the context menu.
- 4 In the context menu, select the item *Properties*.
- 5 In the *Properties* window, select the tab *Drivers*.

In the *Drivers* tab, you can see the version of the installed driver.

6 Operation

6.1 Operation with GetPos

Using the *GetPos* software, you can send position values to your PC, and have the option of saving them in an ASCII text file (.txt) or an Excel-compatible file (.csv) for further processing.

To send position values from the SensorReady 3D to your PC and to process them with *GetPos*, the SensorReady 3D and PC must be correctly connected and configured.

» *Chapter 5.2 Establishing Connection to a PC, page 75*

6.1.1 Installation and Configuration

► Installing *GetPos*

- 1 Create a new folder on your hard drive, to which you would like to install *GetPos*. Name the folder, e.g. *GetPos*.
- 2 Insert the installation CD in the CD drive.
- 3 Open the folder `..\Software\GetPos`.
- 4 Copy the file *GetPos.exe* to the folder on your hard drive which you added in step 1.

GetPos is now installed on your PC.


When *GetPos* is used for the first time, it generates the configuration file *GetPos.cfg*, where the settings you have made are saved. The configuration file is saved in the folder in which the file *GetPos.exe* is stored.

► Identifying COM interface

To configure *GetPos*, you need to identify the COM interface which is used by SensorReady 3D.

- 1 Open the device manager:
 - **Windows 2000, Windows XP:** *Start → Control panel → System → Hardware → Device manager*
 - **Windows Vista:** *Start → right-click on Computer → Device manager*
 - **Windows 7:** *Start → Control panel → System → Device manager*
 - **Windows 8:** Right-click in the left, lower corner of the desktops/start screen → *Device manager*
 - **Windows 10:** Right-click the Windows symbol → *Device manager*

- 2 Expand the column *Ports (COM and LPT)* by double-clicking.

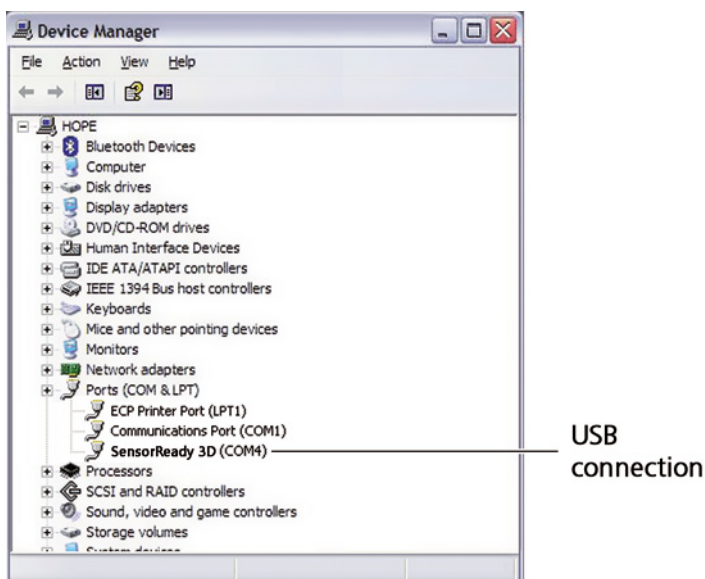The COM interface being used is shown in brackets after the entry *SensorReady 3D*.


The number of the virtual COM interface is retained regardless of which USB port is used on your PC.

► Configuring GetPos

- 1 Start GetPos by double-clicking on the file *GetPos.exe*

When GetPos is started for the first time, the dialogue window *GetPos Setup* appears. Here, you can make all the settings required to establish the connection between SensorReady 3D and GetPos.

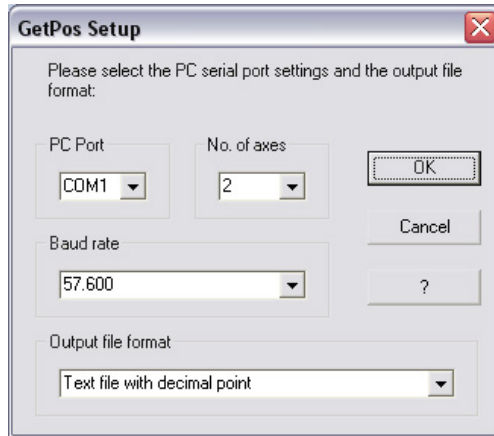Once the settings have been made, the dialogue window *GetPos Setup* no longer appears when GetPos is started. You can, however, adjust all the settings at any time using the *Settings* menu in the menu bar.

» *Chapter 6.1.3 Menu Functions, page 82*

NOTICE

Incorrect data transmission and measurement caused by incorrect configuration.

During configuration of GetPos, always enter the correct interface and the correct number of measuring systems (axes).

- 2 Select the virtual *Port* via which the SensorReady 3D is connected to your PC.


- 3 Select the number of measuring systems (*No. of axes*) which are connected to the SensorReady 3D.
- 4 Select the *File format* to be used when saving your data.
- 5 Confirm your entries by clicking the button *OK*.

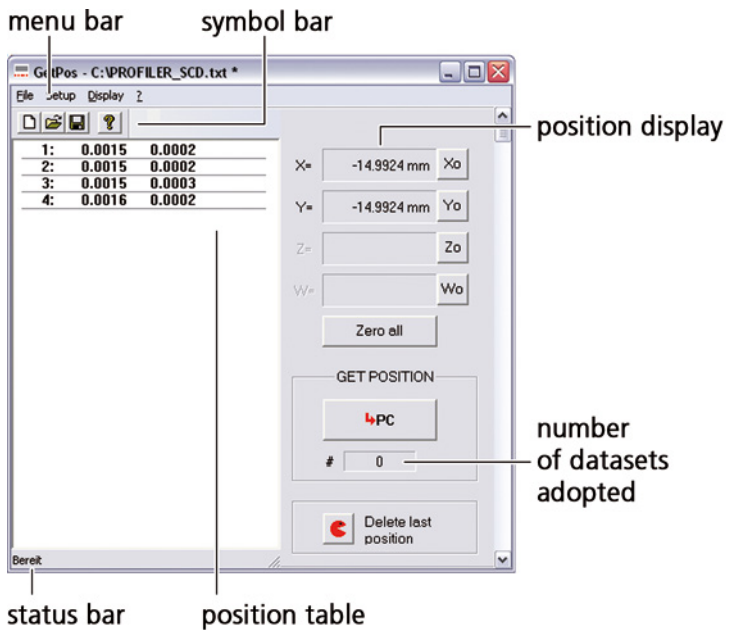GetPos is now configured for communication with the SensorReady 3D.

Selection of the baud rate is only relevant for devices which are connected to the PC via an RS-232 interface. With a USB connection, the baud rate selection has no effect.

6.1.2 User Interface

Button	Function
	Add new position table
	Open saved position table
	Save current position table
	Display version number of GetPos
	Zero X position value
	Zero Y position value
	Zero Z position value
	No function ¹
	Zero position values of all axes
	Adopt current position values in position table
	Delete last adopted position values

¹ Function only available for position displays and measuring interfaces with 4 measuring channels

6.1.3 Menu Functions

Menu	Item	Function
File	New	Add new position table
	Open...	Open saved position table
	Save	Save current position table
	Save as...	Save current position table with a new file name
	Last file	Open last saved position table
	Quit	Quit GetPos
	Settings	COM interface
Baud rate		Select baud rate Note: Selection of the baud rate is only relevant for devices which are connected to the PC via an RS-232 interface. With a USB connection, the baud rate selection has no effect.
File format		Select file format in which the position table is saved
No. of axes		Select the number of measuring systems which are connected to the SensorReady 3D
Language		Select language (German/English)

Menu	Item	Function
Display	Symbol bar	Show/hide symbol bar
	Status bar	Show/hide status bar
?	Info about GetPos...	Display version number of GetPos

6.1.4 Key Layout on PC

Key/key combinations	Function
X Y Z	Zero position value for X, Y, Z
0	Zero position values of all axes
ENTER	Send current position values to PC
CTRL + Backspace	Delete last adopted position value

6.1.5 Sending Position Values to PC

You can send position values to GetPos in different ways:

- Press the button  in GetPos
- Press the *Enter* key on your PC keyboard

When you have sent the position values in one of the ways described, the position values appear in the GetPos position table as a new data set.

6.1.6 Saving Position Tables

You can save position tables in the following file formats:

- Text file with decimal separators (ASCII)
- Text file with comma separators (ASCII)
- Excel-compatible file (.csv)

► Changing file format

When starting GetPos for the first time, specify the different settings, including the file format, in which you will save your position tables. These settings can, however, be adjusted any time using the *Settings* menu in the menu bar.

- 1 In the *Settings* menu in the menu bar, select the preferred file format under the sub-menu *File format*.

The position table is saved in the selected file format the next time it is saved.

► Saving position table

- 1 In the *File* menu in the menu bar, select the item *Save* or click  in the symbol bar.
- 2 Follow the on-screen instructions.

► Saving position table with a new file name

- 1 In the *File* menu, select the item *Save as...*
- 2 Follow the on-screen instructions.

► Opening position table

GetPos automatically adds a new position table with each program start.

- 1 To open a saved position table, click the button  in the symbol bar or select the item *Open...* in the *File* menu in the menu bar.
- 2 Follow the on-screen instructions.

As soon as the position table is visible, you can add additional position values and save the file under a new name or as a different file format.

6.2 Operation via Instruction Set

Using the instruction set, you can read out position values and configure all settings on the SensorReady 3D via the USB interface of the SensorReady 3D.

You will find the SensorReady 3D instruction set in PDF format on the installation CD:

» *SensorReady 3D Instruction Set*

Communication between the SensorReady 3D and the PC takes place using ASCII instructions:

- The entries can be in lower or upper case.
- Instructions are completed with *CR* (carriage return, ASCII: 013_{dec}).
- Return values of the SensorReady 3D are completed with *CR* (carriage return, ASCII: 013_{dec}) as well.

To read out and program the SensorReady 3D using the instruction set, the SensorReady 3D and PC must be correctly connected and configured.

» *Chapter 5.2 Establishing Connection to a PC, page 75*

6.3 Using Position Correction

6.3.1 How Position Correction Works

The SensorReady 3D is equipped with section-by-section linear position correction in order to compensate length errors of your positioning system. The measurements can be corrected using up to 128 reference positions for each axis. The SensorReady 3D calculates position values between the reference positions by linear interpolation.

6.3.2 Configuring Position Correction

Heed the following warnings during configuration of position correction:

NOTICE

Incorrect position correction caused by moving the machine zero or changing the encoder type.

Position correction values are valid as long as the location of the machine zero and the selected encoder type remain unchanged. If, after specifying position correction values, you

- move the machine zero
- mechanically move the reference switches
- approach another reference mark to set the machine zero
- change the function for determining the machine zero or
- change the encoder type

all the position correction values already saved must be deleted and new position correction values must be defined.

Incorrect position correction due to excessive deviation between target and actual value.

The values for *Target Pos* and *Actual Pos* must not deviate from one another by more than 127 μm . Check each value pair before saving it.

6.3.3 Enabling Position Correction

To enable position correction you must approach the machine zero via the reference marking or set it manually once every time you switch the SensorReady 3D on.

7 Cleaning

WARNING

Danger to life and possible device damage due to improper cleaning.

Disconnect the power supply of the SensorReady 3D from the mains before starting cleaning work. Ensure that no liquids enter the interior of the SensorReady 3D or power supply. Protect the connections/ports and control elements of the SensorReady 3D from liquids. Never touch the power supply and socket if your hands are damp or wet.

Danger to life through the spread of substances which are harmful to health.

In the case of possible contamination with corrosive, infectious, radioactive or other substance which is harmful to health, have the SensorReady 3D professionally disinfected/decontaminated.

► **Cleaning SensorReady 3D**

- 1 Disconnect the SensorReady 3D from the mains.
- 2 Clean the SensorReady 3D with a lint-free cloth. With greater contamination, slightly dampen the cleaning cloth with water or a mild solvent-free disinfectant/cleaning agent.
- 3 Allow the SensorReady 3D to dry.
- 4 Connect the SensorReady 3D to the mains.

The SensorReady 3D is now ready to use again.

8 Maintenance

The SensorReady 3D is maintenance-free.

9 Troubleshooting

WARNING

Possible electric shock and equipment damage by connecting external or damaged power supplies.

Only use the 24 V/30 W power supply provided. Immediately replace the power supply if the power supply housing or cable show signs of damage. Only use replacement components from Märzhäuser Sensotech.

Electric shock and equipment damage possible due to unauthorised repair attempts.

Never open the housing of the SensorReady 3D or the power supply provided. Only have repair work carried out by our Customer Service Center. In case of violation, all warranty and liability claims will become void.

The following overview can help you with remedying faults and malfunctions of the SensorReady 3D.

If a fault or malfunction cannot be remedied using the overview, contact our Customer Service Center.

» *Chapter 10 Service, page 92*

Fault	Possible causes
SensorReady 3D cannot be switched on	No power supply
Position values unusual or demonstrably incorrect	Incorrect information about the encoder type or signal period/counting step
	Connection of an incompatible measuring system
	Incorrect cabling
Position correction does not work	Machine zero not set/approached
	Position correction for the affected axis is disabled

Remedy

Check whether the external power supply is properly connected to the port *24V DC/2A max.* on the SensorReady 3D and to the mains.

Check the settings for the encoder type and the signal period/counting step of the measuring system affected in the SensorReady 3D.

» *Chapter 13.2 Encoder Interface, page 98*

Check the specifications of the affected axis.

» *Chapter 4.1 Intended Use, page 67*

- 1 Disconnect the power supply of the SensorReady 3D from the mains.
 - 2 Check whether the measuring system affected is properly connected to the SensorReady 3D via the correct port.
 - 3 Connect the power supply of the SensorReady 3D to the mains.
 - 4 Re-check the position values of the affected axis.
-

When the function *Set machine zero by reference switch* or *Set machine zero by reference mark* has been selected: Approach the reference switches/reference marks of all axes after the SensorReady 3D has been switched on.

Enable the position correction for the affected axis.

10 Service

The staff at our Customer Service Center will be happy to help with support queries, complaints and suggestions about our products and services.

Overview of the services from our Customer Service Center:

- Advice and support
- Support for customer-specific application questions
- Support for putting equipment into service
- Inspection and repair of equipment
- Fault analysis (also online via web conference)
- Delivery of spare parts and replacements
- User training (also online via web conference)

You can reach our Customer Service Center Mondays to Fridays from 08:00 to 16:00.

Tel.: **+49 6441 9116-36**
 Fax: **+49 6441 9116-40**
 E-mail: **service@marzhauser.com**

So that we can help you quickly when you need service, please keep the following information about your SensorReady 3D handy:

Designation	Information on the nameplate
Type	
Article-No.	
Serial-No.	

10.1 Replacement Components

Should you need to replace a component of the SensorReady 3D, please contact our Customer Service Center.

Replacement component	Order no.
MR encoder cable	00-76-650-0804
USB interface cable	00-76-100-2802
External 24 V/30 W power supply incl. 4 socket adapters (EU/GB/USA/AUS)	00-76-203-4801
Operating manual	00-76-901-4830-TD
Installation CD incl. driver, operating manual, software <i>GetPos</i>	00-76-901-4830-CD

11 Taking out of Service

► Taking SensorReady 3D out of service

NOTICE


Possible equipment damage through putting into/taking out of service with existing mains power connection.

Disconnect the SensorReady 3D from the mains before connecting/disconnecting measuring systems, a PC or other peripherals to/from the SensorReady 3D.

- 1 Disconnect the power supply of the SensorReady 3D from the mains.
- 2 Disconnect all connected cables from the SensorReady 3D.
- 3 Clean the SensorReady 3D.
 - » *Chapter 8 Maintenance, page 88*

Taking out of service is now complete.

12 Disposal

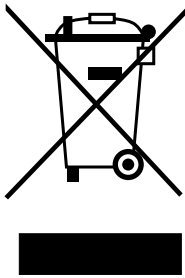The adjacent symbol means that the SensorReady 3D is classified as electrical or electronic equipment in accordance with the EU Directive 2012/19/EU (WEEE) and thus cannot be disposed of with household waste. Since the equipment is used commercially, it may not be disposed of via a communal collection point for electrical and electronic equipment.

Return the SensorReady 3D to Märzhäuser Sensotech for proper disposal free of charge.

All equipment purchased from Märzhäuser Sensotech can be returned to us at any time for free disposal. If you wish to dispose of the equipment yourself, please heed local laws and regulations.

► Disposing of SensorReady 3D

- 1 Take the SensorReady 3D out of service.
» *Chapter 11 Taking out of Service, page 94*
- 2 Pack the SensorReady 3D and accessories properly.
- 3 Send the SensorReady 3D for disposal to:

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Production features

Dimensions	165 x 62 x 29 mm (L x B x H)
Weight	Approx. 400 g (excluding cable)

On request, your SensorReady 3D can also be upgraded for the connection of an incremental measuring system of the type *TTL (non-differential)* to port *Encoder 3*. Please contact our Customer Service Center if necessary:

» *Chapter 10 Service, page 92*

System requirements (when connected to a PC)

Interface	USB interface
Operating system	Windows 2000, Windows XP (32/64 Bit), Windows Vista (32/64 Bit), Windows 7 (32/64 Bit), Windows 8 (32/64 Bit), Windows 10 (32/64 Bit)

Safety functions


Processor monitoring	Reset in case of undervoltage, logic voltage (3.3 V)
Voltage monitoring	Measurement of internal voltages (3.3 V, 5 V, 24 V/DC)
Surge/short-circuit protection (external voltages)	VCC5 + VCCEXT, shutdown is automatically reset after the end of the overload condition

13.2 Encoder Interface

Parameter	MR/ana- logue 5V_{ss}	1V_{ss}	TTL (RS-422)
Position signal	1.5...5 V _{ss} ¹	0.6...1.2 V _{ss} ¹	standard ¹
Reference signal evaluation	–	Yes	Yes
Error signal evalu- ation	On request	On request	On request
Max. input fre- quency	25 MHz	25 MHz	25 MHz ²
Max. interpolation factor	12,868 ³	12,868 ³	4 ⁴
Input resistance	120 Ω	120 Ω	120 Ω
Supply voltage measuring system	5 V ⁵	5 V ⁵	5 V ⁵

¹ Differential transmission

² Pulse interval min. 36 ns

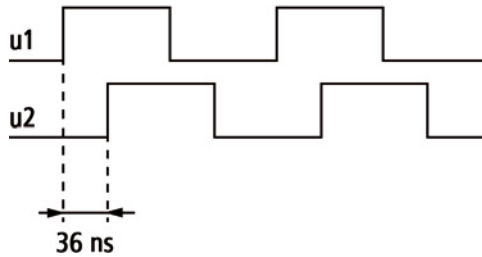³ Arc tangent interpolation

⁴ Quadrature counter

⁵ Max. 500 mA per measuring system, max. 800 mA total

13.3 Ports

At all ports, the cable shield is connected to the housing of the SensorReady 3D.

13.3.1 Encoder 1: HD D-Sub 15 Pin Connector

Pin	Signal	Function
1.2	VCC5	Power supply, 5 V ¹
3	+U1X1	+sine, X axis ²
4	-U1X1	-sine, X axis ²
5	+U2X1	+cosine, X axis ²
6	-U2X1	-cosine, X axis ²
7	REF2	Reference switch input Y axis ³
8	REF1	Reference switch input X axis ³

Pin	Signal	Function
11	+U1X2	+sine, Y axis ²
12	-U1X2	-sine, Y axis ²
13	+U2X2	+cosine, Y axis ²
14	-U2X2	-cosine, Y axis ²
9, 10, 15	GND	Reference potential

¹ Max. 500 mA, total load of all VCC5 ports max. 800 mA

² Output signal: 5V_{ss}, differential transmission, R_L = 120 Ω

³ TTL input, pull-up resistance 1 kΩ → VCC5,
low pass 1 kΩ/10 nF

13.3.2 Encoder 3: D-Sub 9 Pin Connector

Pin assignment for MR/analog 5V_{ss} measuring systems

Pin	Signal	Function
1	-U1	-sine ¹
2	GND	Reference potential
3	-U2	-cosine ¹
4	ERR	Error ²
5	-	-
6	+U1	+sine ¹
7	VCC5	Power supply, 5 V ³
8	+U2	+cosine ¹

Pin	Signal	Function
9	–	–

¹ Output signal: 5V_{ss}, differential transmission, $R_L = 120 \Omega$

² TTL input, pull-down resistance $1 \text{ k}\Omega \rightarrow \text{GND}$,
low pass $1 \text{ k}\Omega/10 \text{ nF}$

³ Max. 500 mA, total load of all VCC5 ports max. 800 mA

Pin assignment for 1V_{ss} measuring systems

Pin	Signal	Function
1	–U1	–sine ¹
2	GND	Reference potential
3	–U2	–cosine ¹
4	ERR	Error ²
5	–U0	–Reference mark measuring system ⁴
6	+U1	+sine ¹
7	VCC5	Power supply, 5 V ³
8	+U2	+cosine ¹
9	+U0	+Reference mark measuring system ⁴

¹ Output signal: 1V_{ss}, differential transmission, $R_L = 120 \Omega$

² TTL input, pull-down resistance $1 \text{ k}\Omega \rightarrow \text{GND}$,
low pass $1 \text{ k}\Omega/10 \text{ nF}$

³ Max. 500 mA, total load of all VCC5 ports max. 800 mA

⁴ Differential transmission, $R_L = 120 \Omega$

Pin assignment for TTL (RS-422) measuring systems

Pin	Signal	Function
1	-U1	Negative U1 ¹
2	GND	Reference potential
3	-U2	Negative U2 ¹
4	ERR	Error ²
5	-U0	-Reference mark measuring system ¹
6	+U1	Positive U1 ¹
7	VCC5	Power supply, 5 V ³
8	+U2	Positive U2 ¹
9	+U0	+Reference mark measuring system ¹

¹ Output signal: TTL, differential transmission, $R_L = 120 \Omega$

² TTL input, pull-down resistance $1 \text{ k}\Omega \rightarrow \text{GND}$,
low pass $1 \text{ k}\Omega/10 \text{ nF}$

³ Max. 500 mA, total load of all VCC5 ports max. 800 mA

13.3.3 USB Interface: 4 Pin Connector, Type B

Pin	Signal	Function
1	VCC1	Power supply (+5 V, max. 500 mA)
2	USBDM	–Data
3	USBDP	+Data
4	GND	Reference potential

¹ internal 100 nF → GND

13.3.4 Feature Connector: HD D-Sub 15 Pin Connector

Pin	Signal	Function
1	–	–
2	–	–
3	GND	Reference potential
4	–	–
5	–	–
6	+IN1	+Reference switch, 0-24 V ¹
7	–IN1	–Reference switch, 0-24 V ¹
8	IN2	Input 1, Clear, 0-24 V1
9	IN3	Input 2, Send, 0-24 V1

Pin	Signal	Function
10	–	–
11	–	Reference potential IN2, IN3
12	NC	Not used
13	–	Reserved for system-internal functions
14	VCC5	Output voltage 5 V ²
15	VCCEXT	Output voltage 24 V ± 10% ³

¹ Logic level: 0–0.8 V = Low; 3.95–28 V = High; secured against reverse polarity, optically isolated

² Max. 500 mA, total load of all VCC5 ports max. 800 mA

³ Max. 0.4 A

13.3.5 24V/DC: 2.1 mm DC Plug

Pin	Signal	Function
Pin, inside	+24 V	Voltage range +8... +28 V, power consumption at +24 V = typ. 110 mA, no load, internal fuse (not replaceable) = 2 A Note: Only use the 24 V/30 W power supply provided.
Housing, external	GND	Shield port (connected to the housing)

13.3.6 Shield: Insert Nut M3

On the insert nut M3 (*shield*), you can set the SensorReady 3D to a specific potential (e.g. protective ground).

The port is not suitable for implementing a standards-compliant PE connection.

13.4 External 24 V/30 W Power Supply (PSAA30R-240)

30W Interchangeable AC Clip Adapter

PSAA30R-XXX-R Series

● 12V/15V/24V/56V ● DoE VI ● Interchangeable AC Clip ● Cable ● Black

The 30W interchangeable AC clip adaptor series is the latest offering from a world leader in charging adapter technology and manufacturing - Phihong Technology. Phihong 30W interchangeable AC clip adaptor series are designed to work with most 12V/15V/24V/56V powered devices and conform to most of the industry standards and meet DoE level VI energy efficiency requirements. We ensure our products meet or exceed all the industry performance and safety requirements.

Features :

Efficiency	DoE Level VI (86.95% avg)
Power Saving	<0.1W @230Vac/50Hz
Protection	SCP, OVP, OCP
Turn On Delay Time	<3s @ 115Vac, max load
Hold-up Time	>10ms @ 115Vac, max load
Rise Time	<100ms @ 115Vac, max load
Isolation	Primary to Secondary 3000Vac, 10mA for 1min.
Insulation	Primary to Secondary >7M ohm 500Vdc
DC Output Cable (12V)	1500±50mm, UL2468, AWG#18
(15V, 24V, 56V)	1500±50mm, UL2468, AWG#24
DC Output Connector	5.5x2.1x10mm, center positive, straight barrel
Dimension	82.7*55*39.1 (mm)
Weight	180 (g)

Electrical Specifications :

Input Voltage	100~240Vac (90~264Vac min/max)
Frequency	50/60Hz Nom (47~63Hz min/max)
Input Current	0.8A (RMS) @ 120Vac max load 0.5A (RMS) @ 240Vac max load
Output Voltage	12V / 15V / 24V / 56V
Output Current	2.5A / 2A / 1.25A / 0.54A
Maximum Output Power	30W continuous
Output Ripple & Noise	≤ 120mV (pk-pk) for 12V output
Line/Load Regulation	±5%
Leakage Current	0.25mA max
Inrush Current	~60A @ 100Vac, 25°C, cold start <100A @ 240Vac, 25°C, cold start

Environmental Specifications :

Operating Temperature	0°C ~ +40°C
Storage Temperature	-20°C ~ +75°C, non-condensing
Operating Humidity	20% ~ 90% RH
Storage Humidity	20% ~ 95% RH, non-condensing

Safety and Compliance :

UL/cUL60950-1, SAAAS/NZS 60950-1:2003, CE approval
FCC Class B, EN55022 Class B, C-Tick approval
ESD : E61000-4-2, Level 4
Radiated Immunity: EN61000-4-3, Level 2
EFT: EN61000-4-4, Level 2
Surge: EN61000-4-5, Level 3
Conducted Immunity: EN61000-4-6, Level 2
Magnetic Immunity: EN61000-4-6, Level 1
Voltage Dips: EN61000-4-11

Note:

- Measurements shall be made with an oscilloscope with 20MHz bandwidth. Output shall be bypassed at a connector with a 0.1uF ceramic capacitor and a 10uF low ESR electrolytic capacitor.
- Efficiency: Measure efficiency at 100%, 75%, 50%, 25% load after burn-in 30 minutes at 100% load, then compute the arithmetic average of these 4 values.
- Interchangeable AC clip: US, Brazil, CN, EU, Korea, UK, India, Argentina, Australia, IEC320 C8

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declares that the following product
 Product name: **SensorReady 3D**
 Article no.: **00-76-901-4830**
 Product type: **USB interface for incremental length measuring systems**

is in conformity with the following European Union harmonisation legislation, when operated in accordance with the operating manual:
EMC directive 2014/30/EU
Low voltage directive 2014/35/EU
RoHS directive 2011/65/EU

The following harmonised standards have been applied:
DIN EN 61010-1:2011-07
DIN EN 61326-1:2013-07

In case of a modification of the product, which was not approved by the manufacturer, this Declaration of Conformity loses its validity.
 This Declaration of Conformity is issued under the sole responsibility of the manufacturer.

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Approved by Uwe Schnitter | 13.09.2016
Issued by Christopher Koch | 13.09.2016 | 10:30

Proved to be QMS conform by Christopher Koch | 13.09.2016
QM document ABS02_QF07_22 | Change index 2 | Page 1 of 1

15 Index

A

Accessories 63

C

Cables 63

Cleaning 88

COM interface 77

Customer Service Center 92

D

Decl. of conformity 65, 108

Delivery contents 63

Disposal 95

Driver 76

E

Encoder cables 63

Encoder interface 98

F

File format 77, 84

G

GetPos 77

Configuration 79

File format 77, 84

Installation 77

Key layout 83

Menu functions 80, 82

Position tables 84

User interface 88

I

Installation requirements 67

Instruction set 85

Intended use 67

M

Maintenance 88

Measuring system

3rd 73

appropriate 67

N

Nameplate 64, 92

O

Operation 77

P

PC 75

Ports 66, 99

Position correction 64, 86

enable/disable 87

Position value

send 81, 83

Power supply 106

Putting into service 71

R

Replacement components 93

S

Safety 67

Safety functions 97

Service 92

Shield 105

Supply voltage 96

Support 92

System requirements 97

T

Taking out of service 94

Technical data 96

Troubleshooting 89

U

User training 92

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An den Fichten 35
35579 Wetzlar (Germany)

Tel.: +49 6441 67106-0
Fax: +49 6441 67106-28

info@marzhauser-st.com
www.marzhauser-st.com

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© 2020 | All rights reserved. Any complete or partial copy requires written authorization from Märzhäuser Sensotech GmbH.
Märzhäuser Sensotech GmbH