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Zubehör 10

Table of Contents

1	Basics of the Operating Manual.....	77
1.1	Symbols and Abbreviations	78
2	Delivery Contents.....	79
2.1	Optional Accessories	80
3	Product Description	81
3.1	Identification.....	81
3.2	Area of Application	81
3.3	Product Structure	82
3.4	Compliance with Directives and Standards	85
4	Safety	86
4.1	Intended Use.....	86
4.2	Installation Requirements	87
4.3	General Safety Notes.....	87
4.4	Warning Note Categories.....	90
4.5	Hazard Symbols Used.....	91
5	Putting Into Service	92
5.1	Setting up Hose Connections	93
5.2	Fitting the Objective Adapter.....	96
5.3	Filling with Immersion Medium	98
5.4	Supplying the Liquid Dispenser with Voltage	100
5.5	Carrying Out a Test Run	101
5.6	Establishing a Connection to a PC	102
6	Configuration	105
6.1	Configuration by Software	105
6.2	Configuration by Instruction Set.....	118
7	Operation	119
7.1	Basic Principles.....	119
7.2	Refilling Immersion Medium.....	120
7.3	Operation in Stand-Alone Mode.....	121
7.4	Operation Using a TANGO Controller	124

7.5	Operation with an External Controller	125
7.6	Operation at a PC (Instruction Set)	129
8	Cleaning	130
8.1	Cleaning the Housing.....	131
8.2	Removing the Immersion Medium.....	131
8.3	Dismantling the Objective Adapter.....	133
9	Maintenance	134
10	Troubleshooting	135
11	Service	138
11.1	Replacement Components	139
12	Taking Out of Service	139
13	Disposal	140
14	Technical Data.....	141
14.1	Power Adapter Recommendation.....	142
15	Index.....	143

1 Basics of the Operating Manual

This operating manual gives you the information which is needed for trouble-free and safe use of the Liquid Dispenser.

You must have fully read and understood this operating manual and especially the safety instructions before you unpack, put into service, operate, clean, take out of service or dispose of the Liquid Dispenser.

This operating manual is directed at people with basic technical knowledge related to handling equipment such as is described here. If you do not have any experience with such equipment, have experienced people help you.


Märzhäuser Wetzlar accepts no liability for damage and operating faults which result from ignoring this operating manual.

1.1 Symbols and Abbreviations

Symbol	Meaning
	Additional information and useful notes
	Start of an instruction
1, 2, 3 ...	Steps within one instruction
»	Cross-reference to another part of this operating manual, another document or a website

Abbreviation	Meaning
Fig.	figure
approx.	approximately, roughly
if nec.	if necessary
incl.	including
LED	light emitting diode
min. / max.	minimum / maximum
s	second
typ.	typically
e.g.	for example

2 Delivery Contents

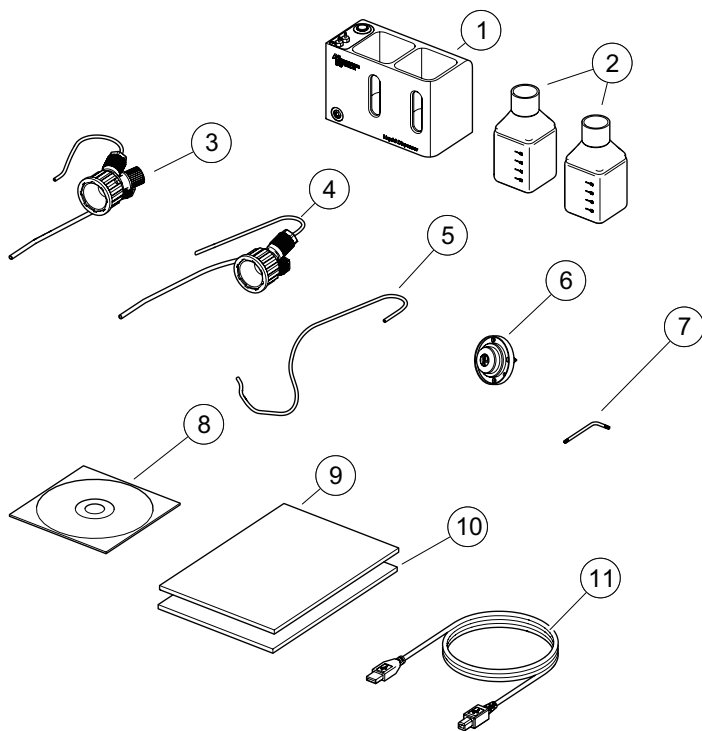
Delivery Contents

No.	Description	Quantity
1	Liquid Dispenser	1
2	Laboratory glass bottle, 100 ml	2
3	Bottle multiple distributor 3x8 incl. PTFE and pneumatic hoses	1
4	Bottle multiple distributor 2x8 incl. PTFE and pneumatic hoses	1
5	PTFE hose	1

No.	Description	Quantity
6	Objective adapter	1
7	TORX wrench T8	1
8	Installation CD	1
9	Operating manual	1
10	Declaration of Conformity	1
11	USB cable	1

2.1 Optional Accessories

Description	Order no.
Wired remote control with start-stop button incl. LED Dimensions: 80 × 50 × 30 mm Cable length: 1.9 m	00-74-300-7000
Connection cable for TANGO controller	00-74-300-6873
Further objective-specific objective adapters	On request
Customer-specific connection cable for external controllers	On request

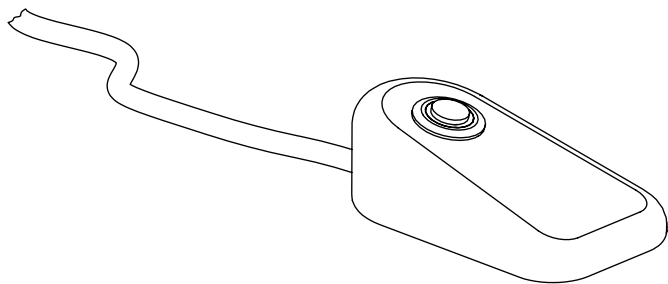Wired remote control

3 Product Description

3.1 Identification

You can find all required information for clear identification of the Liquid Dispenser on the nameplate. The nameplate is attached to the bottom of the Liquid Dispenser.

The objective adapter can be identified on the base of the article number on the bottom.


3.2 Area of Application

The Liquid Dispenser is an automated contamination-free applicator for high-resolution inverted light microscopy. Depending on the version, the Liquid Dispenser has been designed for oil or water immersion. Adaptation to inverted light microscopes from all leading manufacturers is possible using objective-specific objective adapters.

The device can be used both in stand-alone operation and in conjunction with a TANGO controller from Märzhäuser, an external controller or a PC.


Power can be supplied to the Liquid Dispenser either via the USB port or the external I/O interface.

3.3 Product Structure

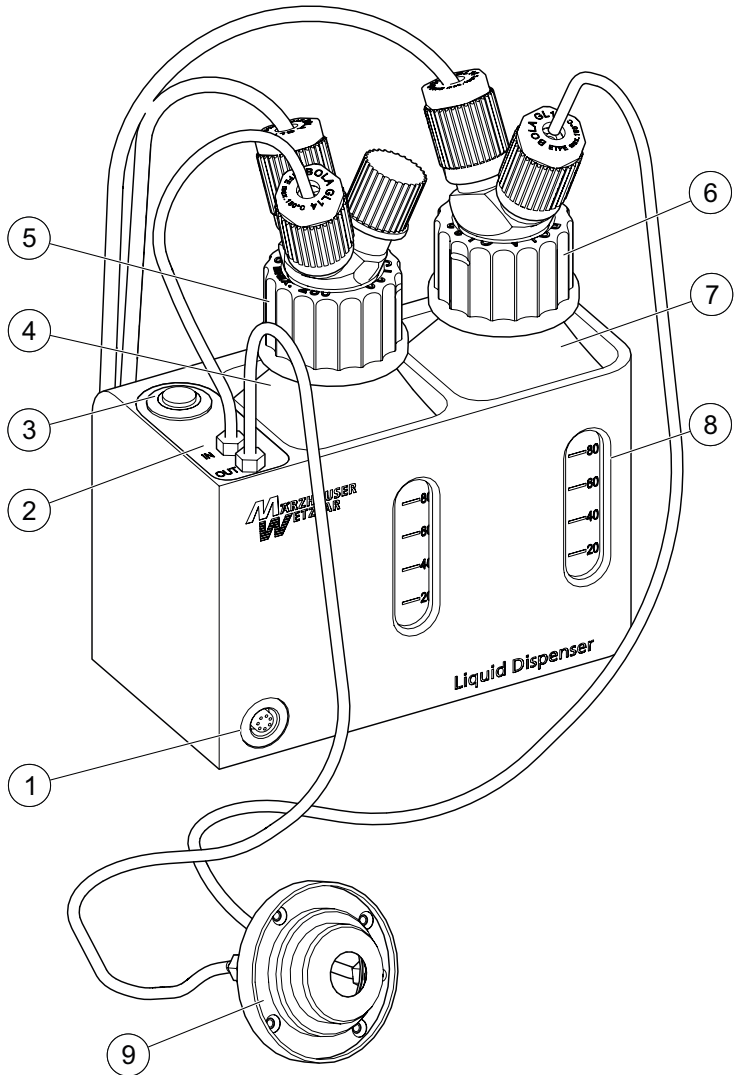Front view

No.	Description	Function
1	7-pin port	Port for remote control
2	Solenoid valve	Controls the dispersion of immersion medium
3	Start/stop button with LED	Operating element and status display
4	Laboratory glass bottle 1	Container for 100 ml unused immersion medium
5	Bottle multiple distributor 3×8	Bottle cap with hose connections
6	Bottle multiple distributor 2×8	Bottle cap with hose connections
7	Laboratory glass bottle 2	Container for 100 ml used immersion medium
8	Viewing window	Filling level display for immersion medium
9	Objective adapter	Objective-specific objective adapter

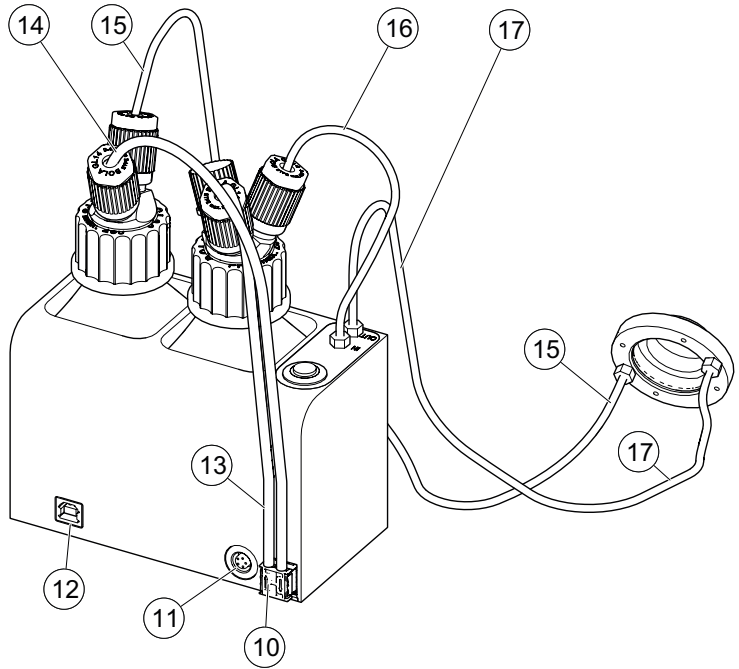Back view

No.	Description	Function
10	Diaphragm pump	Pressure generator for delivering the immersion medium
11	5-pin port	External I/O interface, e.g. connection to the TANGO controller
12	USB port	PC interface (configuration, power supply, control by means of ASCII instruction set)

13	Pneumatic hose (black)	Connection diaphragm pump (↓) to laboratory glass bottle 1
14	Pneumatic hose (black)	Connection diaphragm pump (↓) to laboratory glass bottle 2
15	PTFE hose (white)	Extraction line from objective adapter (<i>OUT</i>) to laboratory glass bottle 2
16	PTFE hose (white)	Delivery line from laboratory glass bottle 1 to the solenoid valve (<i>IN</i>)
17	PTFE hose (white)	Delivery line from the solenoid valve (<i>OUT</i>) to the objective adapter (<i>IN</i>)

3.4 Compliance with Directives and Standards

CE In operation in accordance with the operating manual, the Liquid Dispenser fulfills all the relevant requirements of the following harmonisation regulations of the European Union:

- Machinery Directive 2006/42/EC
- EMC Directive 2014/30/EU
- RoHS Directive 2011/65/EU

The following harmonised standards were applied:

- DIN EN 60204-1:2019
- DIN EN 61326-1:2013
- DIN EN ISO 12100:2011

The complete Declaration of Conformity is included as a separate document with the Liquid Dispenser.

4 Safety

4.1 Intended Use

The Liquid Dispenser

- has been designed and constructed for the contamination-free dispensing of an immersion medium onto the objective of an inverted light microscope.
- depending on the version, is suitable for oil-based immersion media (00-74-300-0000) or for the immersion medium water (00-74-301-0000).
- dispenses fresh immersion medium to the objective adapter and removes used immersion medium again.
- can be operated in manual mode or in interval mode.
- can be operated in stand-alone operation, via the external I/O interface (e.g. with a Märzhäuser TANGO controller), with the aid of a wired remote control or using a PC (configuration software or ASCII instruction set).
- is supplied with power via the USB port (mains adapter or PC) or the external I/O interface (5-pin port) depending on the mode of operation.

The objective adapter is suitable only for fitting to the designated objective.

» *Chapter 3.1 Identification, page 81*

4.2 Installation Requirements

Category	Condition
Installation location	Closed rooms, no explosion hazard, no open electrical ports
Installation height	Max. 2,000 m
Environmental temperature	0...+40 °C
Storage temperature	-20...+70 °C
Cooling	Normal convection
Air humidity	Max. 85 % non-condensing

4.3 General Safety Notes

- Read and make sure you understand the operating manual of the microscope and the peripherals you would like to connect the Liquid Dispenser to.
- Only use the Liquid Dispenser for the purposes which are described in this operating manual.
- It may be necessary to reduce the number of objectives on the objective turret in order to fit the objective adapter.
- It is not possible to pivot the objective turret with the objective adapter fitted. Make sure that the objective turret is not pivoted as long as the objective adapter is fitted.
- Depending on the stage insert or sample holder used, the scan range of your microscope can be restricted with the objective adapter fitted. After you have fitted the objective adapter, check the complete system carefully for possible

collision areas during all operating scenarios (e.g. placement of the stage insert or the sample holder, focussing etc.).

- Always disconnect the equipment you would like to connect to the Liquid Dispenser or disconnect from the Liquid Dispenser from the mains first.
- Always connect the Liquid Dispenser to a power supply which corresponds to the specifications. When using a 5 V USB power adapter, make sure it is compatible with your local mains.
- Set the Liquid Dispenser up in such a way that you can immediately disconnect the power supply without hindrance in case of emergency.
- Always use shielded cables for signal transmission.
- Specific immersion oils can become gummy after a certain time. Find out about the properties of your immersion oil before you use it in the Liquid Dispenser. If necessary, change the immersion medium at regular intervals and clean the Liquid Dispenser and the objective at regular intervals in order to avoid damage. Particular care is required following longer periods of standstill.
- Certain immersion media can contain substances which are harmful to health. For this reason, always heed the safety data sheet provided by the manufacturer or your in-house operating manual and the usual safety precautions for handling chemicals before using your immersion medium.
- The Liquid Dispenser is not protected against corrosive, infectious, toxic, radioactive or other substances which are harmful to health. When handling such substances, ensure that all legal requirements are fulfilled, especially the national regulations on accident prevention.
- Operate the Liquid Dispenser exclusively with the immersion medium it has been designed for.
- Ensure that no chips, liquids, objects or other substances enter the interior of the Liquid Dispenser.

- Protect the Liquid Dispenser from dust and dirt as far as possible.
- Ensure that the Liquid Dispenser, in conjunction with your application, meets the applicable safety regulations and legal requirements.
- Disconnect all cables from the Liquid Dispenser immediately if you find any malfunctions, faults, safety defects or leaks on the objective adapter or at hose connections. In this case, secure the Liquid Dispenser against being put back into operation.
- If the Liquid Dispenser hoses or cables connected are damaged, replace these without delay.
- Never open the housing of the Liquid Dispenser. In case of violation, all warranty and liability claims will become void.
- Only have repair work carried out by our Customer Service Center.
 - » *Chapter 11 Service, page 138*
- If you lend or sell the Liquid Dispenser you must include all the documents provided.

4.4 Warning Note Categories

Warning notes are divided into the following categories:

WARNING

Hazardous situation that can cause serious injury or death if not avoided

CAUTION

Hazardous situation that can cause light to moderate injury if not avoided

NOTICE

Warning of possible damage or malfunction

4.5 Hazard Symbols Used

The following hazard symbols are used on the Liquid Dispenser:

Symbol	Meaning
	Warning of a general hazard which can be avoided by observing the operating manual

The following hazard symbols are used within warning notes:

Symbol	Meaning
	Warning of fire hazard
	Warning of hazards caused by substances which are harmful to health
	Warning of cutting injuries

5 Putting Into Service

Read and make sure you understand the *Safety* chapter before you start set-up and putting into service.

» *Chapter 4 Safety, page 86*

WARNING

Possible health risk and equipment damage by putting into service too early.


If you use the Liquid Dispenser in regular operation before you have completed putting into service properly, this can result in health risks and equipment damage.

- ▶ Carry out all the steps related to putting into service described in the chapters 5.1 to 5.5.
-


5.1 Setting up Hose Connections

► Setting up hose connections


- 1 Take the lids off the laboratory glass bottles.
- 2 Place the laboratory glass bottles in the Liquid Dispenser and make sure that the scales of the laboratory glass bottles are pointing towards the Liquid Dispenser viewing windows.

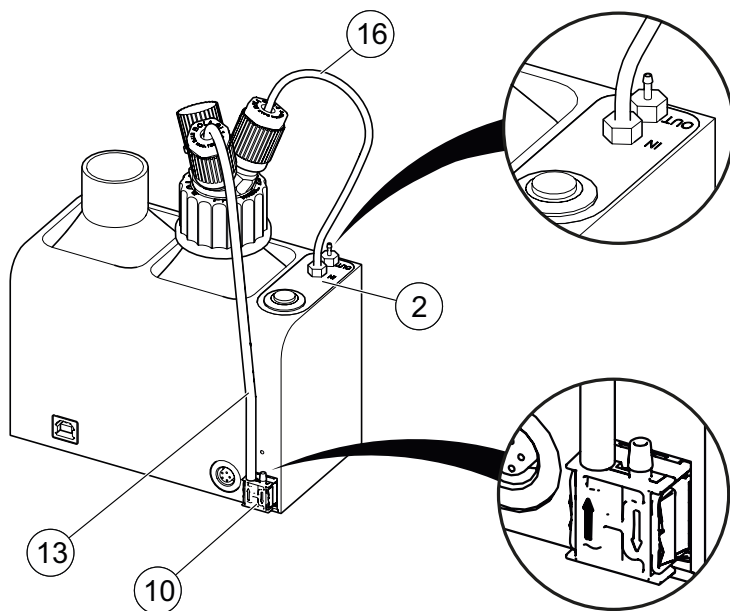- 1 Fix the bottle multiple distributor 3×8 to laboratory glass bottle 1. To do this, hold the hose connections tight with one hand while turning the bottle multiple distributor clockwise with the other hand.
- 2 Fit the PTFE hose ⑯ to the *IN* connection of the solenoid valve ②.
- 3 Fit the pneumatic hose ⑬ to the left connection ↑ of the membrane pump ⑩.

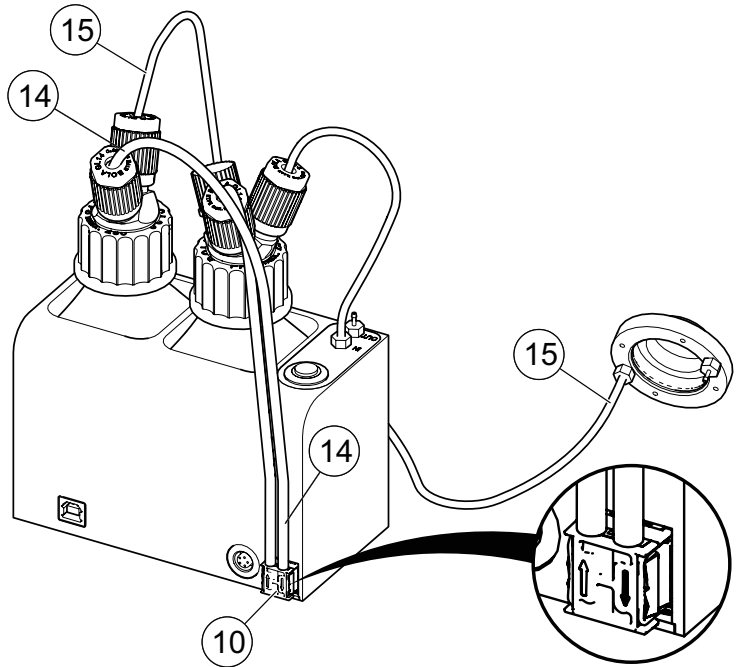- 4 Fix the bottle multiple distributor 2×8 to laboratory glass bottle 2. To do this, hold the hose connections tight with one hand while turning the bottle multiple distributor clockwise with the other hand.
- 5 Fit the pneumatic hose ⑭ to the right connection ↑ of the membrane pump ⑩.
- 6 Fit the PTFE hose ⑮ to the *OUT* connection of the objective adapter.

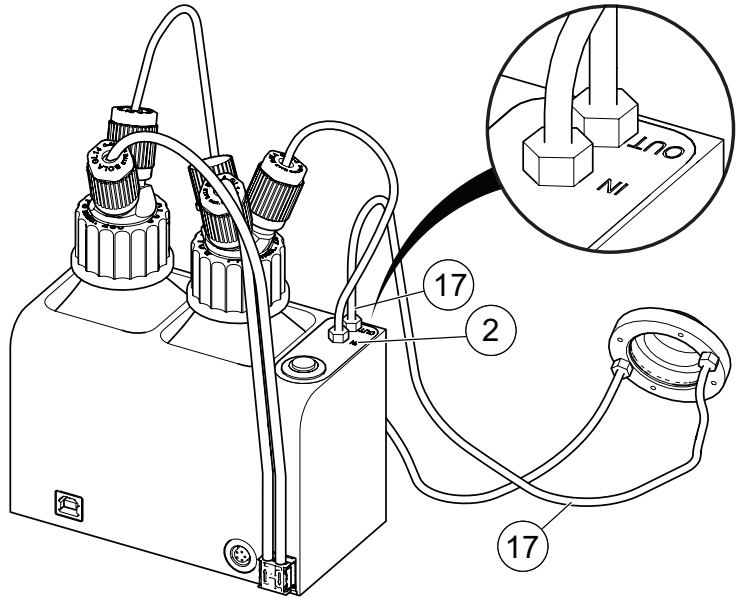- 1 Fit the PTFE hose ⑰ to the *OUT* connection of the solenoid valve ②.
- 2 Fit the other end of the PTFE hose ⑰ to the *IN* connection of the objective adapter.

The hose connections have been set up. Now fit the objective adapter.

5.2 Fitting the Objective Adapter

NOTICE

Possible device damage due to incorrect fitting.

Leaks can occur if an unsuitable objective adapter is used or if the adapter is not fitted correctly. In this case, the immersion medium can leak and cause damage to the microscope.

- ▶ Always use the suitable objective adapter when working with different objectives.
 - » *Chapter 3.1 Identification, page 81*
- ▶ Carry out a test run before starting standard operation.
 - » *Chapter 5.5 Carrying Out a Test Run, page 101*

Possible device damage due to incorrect operation.

The objective turret can no longer be pivoted after the objective adapter has been fitted and can lead to collisions.

- ▶ Make sure that the objective turret is no longer pivoted after the objective adapter has been fitted.

▶ Fitting the objective adapter to the objective

- 1 Remove neighbouring objectives if necessary to be able to fit the objective adapter to the required objective without any obstructions.
- 2 Fit the objective adapter to your microscope objective.


Note: Depending on the working distance, make sure that the objective adapter is approx. 0.05–0.2 mm above the top edge of the objective and aligned horizontally.

The screws of the objective adapter are not tightened when the equipment is delivered so that the objective adapter can be fitted easily on the objective. The screws must be tightened after fitting.


If immersion medium leaks at the objective adapter screws

during the test run or standard operation, the screws must be tightened again.

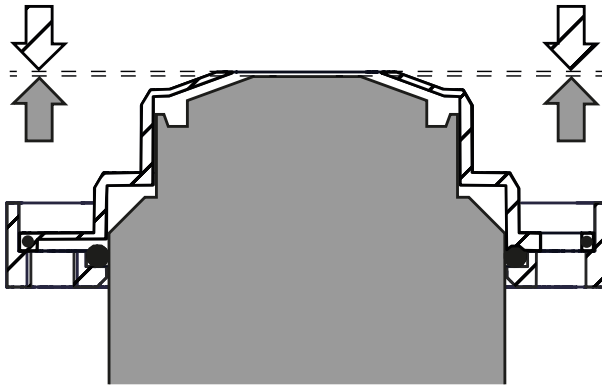 Objective adapter  Objective  Sealing ring

- 3 Use the Torx key T8 to tighten the screws ① to ④ evenly and crosswise. This clamps the objective adapter tight to the objective.

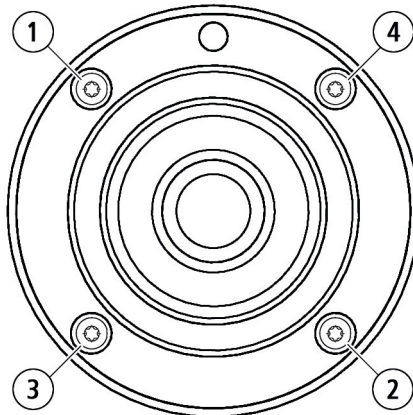The objective adapter has been fitted.

- 4 Check the complete system carefully for possible collision areas during all operating scenarios (e.g. placement of the stage insert or the sample holder, focussing etc.).

If the use of the objective adapter does not cause a collision during operation, you can now fill the immersion medium into the Liquid Dispenser.

» *Chapter 5.3 Filling with Immersion Medium, page 98*

5.3 Filling with Immersion Medium

CAUTION

Danger of cutting by damaged laboratory glass bottles.

The laboratory bottles are made of glass and can be dropped and broken if not handled carefully. In this case, the immersion medium can leak.

- ▶ Make sure you handle the laboratory glass bottles carefully.

NOTICE

Equipment damage by using the wrong immersion medium.

Depending on the version, your Liquid Dispenser has been designed for oil-based immersion media or for the immersion medium water. Use of the wrong immersion medium can make the Liquid Dispenser useless for use of the actual immersion medium.

- ▶ Operate the Liquid Dispenser exclusively with the immersion medium it has been designed for.

» *Chapter 4.1 Intended Use, page 86*

- ▶ Note the specifications of the objective manufacturer.
- ▶ Always use unused immersion medium.

Device damage caused by the immersion medium becoming gummy.

Specific immersion oils can become gummy after a certain time. Particular care is required following longer periods of standstill.

- ▶ Find out about the properties of your immersion oil before you use it in the Liquid Dispenser.
 - ▶ If necessary, change the immersion medium at regular intervals and clean the Liquid Dispenser and the objective at regular intervals in order to avoid damage.
-

▶ Filling with Immersion Medium

- 1 Disconnect the Liquid Dispenser from the power supply.
- 2 Remove the bottle multiple distributor from laboratory glass bottle 1. To do this, hold the hose connections tight with one hand while turning the bottle multiple distributor counterclockwise with the other hand.
- 3 Fill max. 100 ml immersion medium into laboratory glass bottle 1 Use a funnel if necessary.
- 4 Fix the bottle multiple distributor to laboratory glass bottle 1. To do this, hold the hose connections tight with one hand while turning the bottle multiple distributor clockwise with the other hand.

The immersion medium has been filled. The Liquid Dispenser can now be connected to the power supply.

5.4 Supplying the Liquid Dispenser with Voltage

WARNING

Fire hazard and possible equipment damage caused by connection to a non-specified power supply.

The Liquid Dispenser has been designed for a power supply via USB or the 5-pin port. If the specifications are not heeded, fire hazard and equipment damage can be the result.

- ▶ Make sure that your power supply complies with the specifications of the Liquid Dispenser.
 - » *Chapter 14 Technical Data, page 141*

▶ **Supplying the Liquid Dispenser with voltage**

- 1 Connect the Liquid Dispenser to a 5 V power supply. Power can be supplied to the Liquid Dispenser in different ways:
 - Connection via the USB cable included to a PC
 - Connection via an (optional) cable to a Märzhäuser TANGO controller or an external controller
 - Connection via the USB cable included to a USB power adapter

The LED on the start/stop button flashes.

» *Chapter 7.1 Basic Principles, page 119*

The Liquid Dispenser is connected to a power supply. Now carry out a test run.

Excess voltage is not possible when the USB port and a correctly connected controller are used at the same time.

5.5 Carrying Out a Test Run

NOTICE

Possible device damage caused by immersion oil leaking.

Incorrect fitting can cause leaks from the objective adapter. In this case, the immersion medium can leak and cause damage to the microscope.

- ▶ Make sure that no immersion medium can drip down the objective during the test run.

▶ Carrying Out a Test Run

- 1 Press the start/stop button to start dispensing of the immersion medium.
- 2 Wait until the LED goes off.
- 3 Repeat the procedure until the objective adapter has been flooded with medium.
- 4 Make sure that
 - the hose connections are connected properly and
 - that no immersion medium is leaking on the bottom of the objective adapter.
- 5 Let the Liquid Dispenser dispense immersion medium until some immersion medium has been collected in laboratory glass bottle 2.
- 6 Make sure that
 - the hose connections are connected properly,
 - no immersion medium is leaking on the bottom of the objective adapter and
 - that the immersion medium is being extracted properly.

Putting into service of the Liquid Dispenser is now complete. The Liquid Dispenser is now ready for operation.

5.6 Establishing a Connection to a PC

To configure the Liquid Dispenser, connection to a PC and installation of a driver are necessary.

In stand-alone operation you can only change the dispensing duration and the operating mode of the Liquid Dispenser.

» *Chapter 7.3 Operation in Stand-Alone Mode, page 121*

5.6.1 PC System Requirements

Hardware

PC Pentium 4 and higher, CD drive, 1 free USB port (at least USB 2.0), 15 MB free hard drive capacity, monitor (min. resolution 1024 × 768 pixels)

Operating systems supported

All Windows versions Windows 2000 and higher

5.6.2 Installing the Driver

► Installing the driver

Driver installation is explained using Windows 7 as an example. If you are using Windows Embedded, make sure that all required drivers are installed and especially that the *usbser.sys* driver file is in place.

Observe the following notes, regardless of the version of your operating system:

- Do not connect to *Windows Update*.
- Continue installing the driver despite the message about the missing *Windows-Logo-Test*.
- If necessary, agree to the change query of the user account control.

- 1 Switch the PC on.
- 2 Insert the installation CD in the CD drive.
- 3 Connect the Liquid Dispenser to a USB port on your PC.

Windows recognises the Liquid Dispenser as a new device, but the driver is not installed automatically.

- 4 Open the device manager and expand the category *Other devices*.
- 5 Right-click the entry *Liquid Dispenser*.
- 6 In the context menu, select the item *Update driver software*.
- 7 Select the button *Browse for driver software on your computer*.
- 8 Click the button *Browse...* and navigate on the CD to
 - `\Driver Package\Driver_Liquid_Dispenser_x86`
(32 bit operating system) or
 - `\Driver Package\Driver_Liquid_Dispenser_x64`
(64 bit operating system)
- 9 Click the button *Next* to install the driver.

The device manager signals the successful installation and shows the new (virtual) COM port in the title bar of the window.

Driver installation is now complete.

The Liquid Dispenser can now be configured and activated by instruction set.

To configure the Liquid Dispenser using the convenient user interface, install the *Liquid Dispenser Configuration Tool*.

» *Chapter 5.6.3 Installing the Liquid Dispenser Configuration Tool, page 104*

5.6.3 Installing the Liquid Dispenser Configuration Tool

► Installing the Liquid Dispenser Configuration Tool

- 1 Create a new folder on your hard drive, to which you would like to install the *Liquid Dispenser Configuration Tool*. Name the folder e.g. *Liquid Dispenser*.
- 2 Insert the installation CD in the CD drive.
- 3 Open the folder *Configuration Tool* on the installation CD.
- 4 Copy the file *Liquid_Dispenser_Configuration_Tool_X.XX.exe* to the folder you created in step 1.

The *Liquid Dispenser Configuration Tool* is now installed on your PC.

6 Configuration

To configure the Liquid Dispenser, connection to a PC is necessary.

» *Chapter 5.6 Establishing a Connection to a PC, page 102*

The dispensing duration and operating mode can also be configured in stand-alone operation without a PC.

» *Chapter 7.3 Operation in Stand-Alone Mode, page 121*


6.1 Configuration by Software

All Liquid Dispenser settings can be made with the aid of the *Liquid Dispenser Configuration Tool*.

- 1 Install the *Liquid Dispenser Configuration Tool* on your PC.
» *Chapter 5.6.3 Installing the Liquid Dispenser Configuration Tool, page 104*
- 2 Connect the Liquid Dispenser to the PC.
» *Chapter 5.6 Establishing a Connection to a PC, page 102*
- 3 Start the *Liquid Dispenser Configuration Tool* by double-clicking *Liquid_Dispenser_Configuration_Tool_X.XX.exe*.

All the functions of the *Liquid Dispenser Configuration Tool* are described in the following chapters.

6.1.1 Overview – Liquid Dispenser Configuration Tool

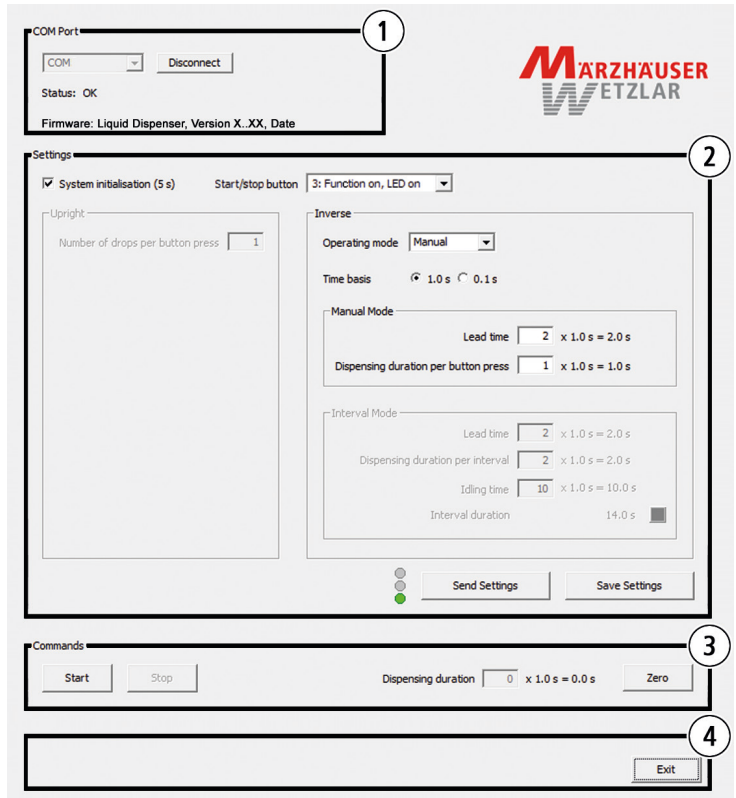The *Liquid Dispenser Configuration Tool* is divided into the following 4 areas:

- ① COM Port
- ② Settings
- ③ Commands
- ④ Quit Program

6.1.2 Area 1 – COM Port

Description	Function
Drop-down menu <i>COM</i>	Select the COM port that the Liquid Dispenser is connected to
Disconnect	Disconnects an existing connection
Status	Displays the connection status (<i>OK</i> or <i>Error</i> with number) Note: All error numbers are described in the Instruction set of the Liquid Dispenser. » <i>Instruction set Liquid Dispenser</i>
Firmware	Version of the firmware installed on the Liquid Dispenser

As soon as the *Liquid Dispenser Configuration Tool* has been connected with the Liquid Dispenser the current device settings are loaded from the Liquid Dispenser and the firmware version is displayed.

The settings are saved directly in the Liquid Dispenser. No data are saved on the PC.

6.1.3 Area 2 – Settings

Manual operating mode

Description

Function

System initialisation (5 s)


Enable/disable system initialisation
 If system initialisation is enabled, the diaphragm pump is switched on immediately for 5 s after the Liquid Dispenser has been switched on. This creates the pressure required to dispense the immersion medium.
 If system initialisation is disabled, the initial dispensing of the immersion medium is delayed.

Start/stop button (0 to 3)	<p>Configure start/stop button</p> <p>Depending on the option selected, the functionality and LED of the start/stop button can be activated or deactivated.</p> <p>Function on: The function of the start/stop button is activated. The start/stop button can be used to operate the Liquid Dispenser.</p> <p>Function off: The function of the start/stop button is deactivated. The Liquid Dispenser can only be operated via controller/PC.</p> <p>LED on: The LED of the start/stop button is activated and indicates the current operating mode of the Liquid Dispenser.</p> <p>LED off: The LED of the start/stop button is deactivated.</p>
Operating mode	<p>Select operating mode</p> <p><i>Manual:</i> When the start/stop button is pressed, immersion medium is dispensed depending on the set <i>Time base</i> and the <i>Dispensing duration per button press</i>.</p> <p>Note: The external I/O interface can also be used in this mode.</p>
Time base (1.0 s or 0.1 s)	<p>Select time base</p> <p>The time base serves as a multiplier for all the following values.</p>

Lead time (0 to 600)	<p>Enter the lead time</p> <p>By means of the lead time, pressure can be built up in the laboratory glass bottle prior to every dispensing of immersion medium. During the lead time only the diaphragm pump is switched on, the solenoid valve remains closed. This ensures that there is always sufficient pressure in the laboratory glass bottle and that immersion medium can be dispensed of safely from the beginning.</p> <p>Duration of pressure build-up = <i>Lead time</i> × <i>Time base</i></p>
Dispensing duration per button press (1 to 6,000)	<p>Enter the dispensing duration per button press</p> <p>The values <i>Dispensing duration per button press</i> and <i>Time base</i> are used to define how long immersion medium is dispensed when the start/stop button is pressed.</p> <p>Dispensing duration = <i>Dispensing duration per button press</i> × <i>Time base</i></p>
Status display (signal light)	<p>Red: The current settings have not been transmitted to the Liquid Dispenser yet.</p> <p>Yellow: The current settings have been transmitted and saved temporarily in the Liquid Dispenser.</p> <p>Green: The current settings correspond to the settings saved in the Liquid Dispenser.</p>
Send settings	<p>Transmit the current settings temporarily to the Liquid Dispenser (no permanent storage in the device)</p>

Save settings	Transmit the current settings to the Liquid Dispenser and save them permanently in the device
---------------	---

Interval operating mode

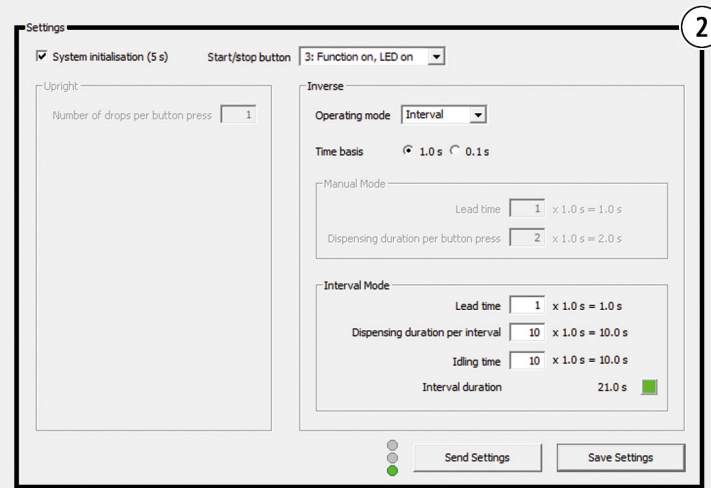Description

Function

System initialisation (5 s)

Enable/disable system initialisation
If system initialisation is enabled, the diaphragm pump is switched on immediately for 5 s after the Liquid Dispenser has been switched on. This creates the pressure required to dispense the immersion medium.

If system initialisation is disabled, the initial dispensing of the immersion medium is delayed.


Start/stop button (0 to 3)	<p>Configure start/stop button</p> <p>Depending on the option selected, the functionality and LED of the start/stop button can be activated or deactivated.</p> <p>Function on: The function of the start/stop button is activated. The start/stop button can be used to operate the Liquid Dispenser.</p> <p>Function off: The function of the start/stop button is deactivated. The Liquid Dispenser can only be operated via controller/PC.</p> <p>LED on: The LED of the start/stop button is activated and indicates the current operating mode of the Liquid Dispenser.</p> <p>LED off: The LED of the start/stop button is deactivated.</p>
Operating mode	<p>Select operating mode</p> <p><i>Interval:</i> When the start/stop button is pressed, immersion medium is dispensed continually depending on the values <i>Time base</i>, <i>Lead time</i>, <i>Dispensing duration per interval</i> and <i>Idling time</i>.</p> <p>Note: In interval mode, the external I/O interface is inactive and thus cannot be used.</p>
Time base (1.0 s or 0.1 s)	<p>Select time base</p> <p>The time base serves as a multiplier for all the following values.</p>

Lead time	<p>Enter the lead time</p> <p>By means of the lead time, pressure can be built up in the laboratory glass bottle prior to every dispensing of immersion medium. During the lead time only the diaphragm pump is switched on, the solenoid valve remains closed. This ensures that there is always sufficient pressure in the laboratory glass bottle and that immersion medium can be dispensed of safely from the beginning.</p> <p>Duration of pressure build-up = <i>Lead time</i> × <i>Time base</i></p>
Dispensing time per interval	<p>Enter the dispensing duration per interval</p> <p>The values <i>Dispensing duration per interval</i> and <i>Time basis</i> are used to define how long immersion medium is dispensed per interval.</p> <p>Dispensing duration immersion medium = Dispensing duration per interval × Time basis</p>
Idle time	<p>Enter the idling time</p> <p>The idling time is the time until the next interval begins.</p>

Interval time

The interval duration is the time which the immersion medium takes place automatically. Immersion medium is dispensed continually.

Time pattern shown on the base of an interval duration of 60 s:
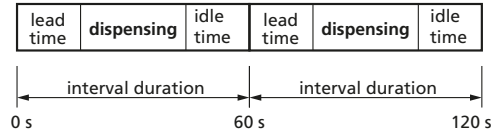
$$\text{Interval duration} = (\text{Lead time} + \text{Dispensing duration per interval} + \text{Idle time}) \times \text{Time base}$$

Status display interval duration (square)

Red: The calculated interval duration is outside the valid value range.


Green: The calculated interval duration is within the valid value range.

Valid value ranges:

Time base 1.0 s: 2.0 s to 6,000 s

Time base 0.1 s: 0.2 s to 600 s


6.1.4 Area 3 – Commands

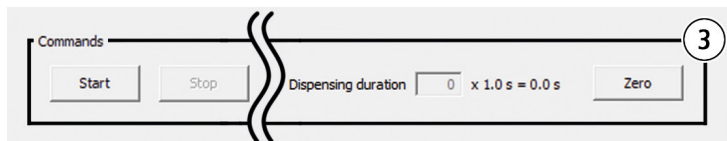With the aid of the functions in the *Commands* area, the Liquid Dispenser can be activated by the software.

Description	Function
Start	<p>Start dispensing of immersion medium using the current settings</p> <p>Note: This function is only available if the current settings have been transmitted to the Liquid Dispenser or saved in the Liquid Dispenser.</p>
Stop	<p>Stop dispensing of the immersion medium immediately</p> <p>The diaphragm pump is switched off and the solenoid valve closed. All processes to be executed are ended.</p> <p>Note: This function is only available if immersion medium is currently being dispensed or interval dispensing is enabled.</p>
Dispensing duration	<p>This counter defines how long immersion medium has already been dispensed, depending on the <i>Time base</i> set.</p> <p>Dispensing immersion medium = $Dispensing\ duration \times Time\ base$</p>
Zero	Reset counter status <i>Dispensing duration</i>

6.1.5 Area 4 – Quit Program

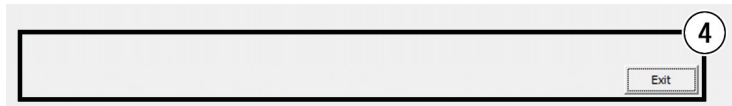
Description	Function
Exit	Quit program If a connection exists via the COM interface, this is disconnected automatically when the program is quit.

6.1.6 Identifying COM Port

This is how to identify the COM port being used by the Liquid Dispenser.

- Open the device manager:
 - Windows 2000, Windows XP:** *Start → Control panel → System → Hardware → Device manager*
 - Windows Vista:** *Start → right-click on Computer → Device manager*
 - Windows 7:** *Start → Control panel → System → Device manager*
 - Windows 8:** Right-click in the left, lower corner of the desktop → *Device manager*
 - Windows 8.1:** Move the mouse pointer to the bottom right-hand corner of the desktop → enter *Device manager* in the search box → *Device manager*
 - Windows 10:** Right-click the Windows symbol → *Device manager*

- Expand the category *Ports (COM & LPT)*.

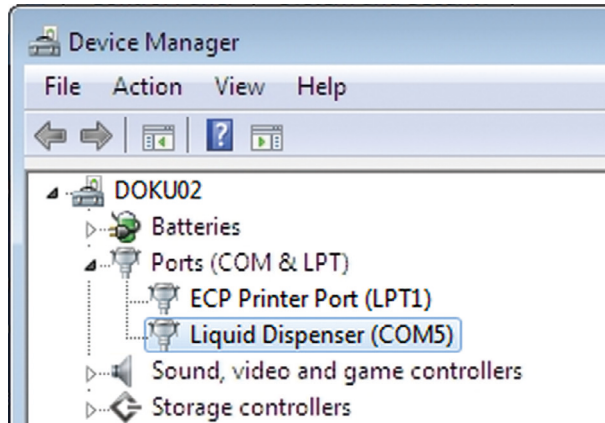The COM port being used is shown in brackets after the entry *Liquid Dispenser*.

The number of the virtual COM port is retained regardless of which USB port is used on your PC.

6.2 Configuration by Instruction Set

All the Liquid Dispenser settings can be configured and immersion medium dispensed by instruction set (ASCII instructions) via the USB interface.

You will find the instruction set on the installation CD:
 » *Instruction Set Liquid Dispenser*

7 Operation

7.1 Basic Principles

Factory settings

The Liquid Dispenser is delivered with the following factory settings:

Parameter	Value	Unit
System initialisation (5 s)	Disabled	–
Time base	1	Seconds
Lead time	0	Factor
Dispensing duration	5	Factor
Idle time	55	Factor
Operating mode	Manual	–

To change to the factory settings, configure the Liquid Dispenser with the aid of a PC or the start/stop button.
 » *Chapter 6 Configuration, page 105*

Switching the Liquid Dispenser on

The Liquid Dispenser is switched on as soon as it is connected to a power supply.

Identifying the operating mode

As soon as the Liquid Dispenser is connected to the power supply, the LED on the start/stop button indicates the current operating mode.

LED	Operating mode
Flashes 2×	Manual mode
Flashes 3×	Interval mode

7.2 Refilling Immersion Medium

- 1 Stop operation of the Liquid Dispenser.
- 2 Disconnect the Liquid Dispenser from the power supply.
- 3 Remove the bottle multiple distributors from the laboratory glass bottles 1 and 2 as follows: Hold the hose connections tight with one hand while turning the bottle multiple distributor counterclockwise with the other hand.
- 4 Dispose of the immersion medium from laboratory glass bottle 2 in accordance with manufacturer's instructions.
- 5 Fill laboratory glass bottle 1 up to the maximum filling height with new immersion medium of the type you were using before.
- 6 Fix the bottle multiple distributors to laboratory glass bottles 1 and 2 as follows: Hold the hose connections tight with one hand while turning the bottle multiple distributor clockwise with the other hand.
- 7 Now connect the Liquid Dispenser to the power supply again.

The Liquid Dispenser is now ready for operation again.

7.3 Operation in Stand-Alone Mode

To operate the Liquid Dispenser in stand-alone mode, it simply needs to be connected to the power supply. The value *Dispensing duration per button press* can be changed using the start/stop button and you can switch between manual mode and interval mode. The *Time base* and the *Lead time* can only be configured using the *Liquid Dispenser Configuration Tool* or the instruction set.

The wired remote control can also be used instead of the start/stop button.

» Chapter 2.1 *Optional Accessories*, page 80

7.3.1 Manual Mode

► Starting/stopping the dispensing of immersion medium

Function	Start/stop button	LED
Start dispensing	Press 1x	Lights up
Lead time	–	Goes out
Dispensing immersion medium	–	Is lit constantly
Interrupt the process if necessary	Press 1x	Goes out
Dispensing completed	–	Goes out

► **Changing the dispensing duration of the immersion medium**

Function	Start/stop button	LED
Start configuration mode	Keep pressed for approx. 5 s	Lights up constantly at first, then flashes slowly
–	Release	Goes out
Set dispensing duration	Press 1× per time unit ¹ (e.g.: press 5× = 5 time units)	Lights up every time button is pressed
Save entry	Keep pressed for approx. 5 s	Lights up constantly at first, then goes out
End configuration mode	Release	–

¹ 1.0 s (factory setting) or 0.1 s

7.3.2 Interval Mode

The interval mode of the Liquid Dispenser makes it possible to dispense the immersion medium at regular intervals.

If interval mode is enabled, this is indicated by the LED flashing continuously. The LED flashes at an interval of approx. 3 seconds (1.5 s LED on, 1.5 s LED off).

► Starting/stopping interval mode

Function	Start/stop button	LED
Start interval mode	Press 1×	Lights up
Lead time	–	Flashes slowly
Dispensing immersion medium	–	Is lit constantly
Idling to end of interval	–	Flashes slowly
The set interval is repeated until interval mode is stopped.		
Stop interval mode	Press 1×	Goes out

► Changing interval mode settings

In order to change the interval mode settings, use the *Liquid Dispenser Configuration Tool* or the instruction set.

- » Chapter 6.1 Configuration by Software, page 105
- » Chapter 6.2 Configuration by Instruction Set, page 118

7.3.3 Changing the Operating Mode

Function	Start/stop button	LED
Change the operating mode	Keep pressed for approx. 20 s	Lights up, flashes slowly, flashes quickly
End mode change	Release	Goes out

7.4 Operation Using a TANGO Controller

► Connecting the Liquid Dispenser with the TANGO controller

- 1 Disconnect the TANGO controller and the Liquid Dispenser from the power supply.
- 2 Connect the Liquid Dispenser and the TANGO controller using the optionally available AUX I/O cable.
- 3 Switch the TANGO controller on.

You can now operate the Liquid Dispenser using the TANGO instruction set from firmware 1.66 and higher.

» *TANGO Instruction Set*

7.5 Operation with an External Controller

The Liquid Dispenser can be operated with an external controller via the 5-pin port and be activated with the aid of various signals.

7.5.1 Signal Characteristics

Signal	Voltage	Polarity	Pulse width	Frequency
Input DROP_IN	3.3 V ...5 V	High active	min. 2.5 ms	max. 100 Hz
Input PUMP_IN	3.3 V ...5 V	High active	min. 100 ms	–
Output FEEDBACK_ OUT	5 V	High active	2 ms	–


The time base can be configured using the Liquid Dispenser *Configuration Tool* or the instruction set.

- » Chapter 6.1 *Configuration by Software*, page 105
- » Chapter 6.2 *Configuration by Instruction Set*, page 118

7.5.2 Signal Description

Signal at the input **DROP_IN**

By means of the signal *DROP_IN*, immersion medium can be dispensed depending on the set time base 0.1 s or 1.0 s. For every impulse received (signal = "High"), immersion medium is dispensed for 0.1 s or 1.0 s. Several successive impulses are added up and processed in sequence.

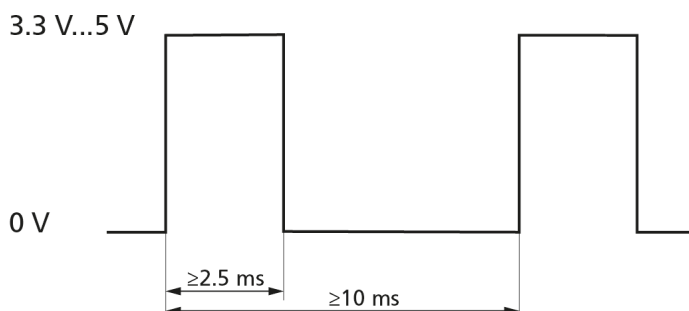
Signal	Action
High	Dispersion of immersion medium for 0.1 s or 1.0 s
Low	–

In order to dispense immersion medium for 5 s with the time base 1 s selected, for example, 5 impulses must be transmitted in succession to the Liquid Dispenser. During impulse transmission the times shown in the diagram must be kept.

Signal at the input PUMP_IN


The *PUMP_IN* signal offers the possibility of switching the diaphragm pump of the Liquid Dispenser on and off manually. Manual switch-on of the diaphragm pump can guarantee that there is always sufficient pressure in the laboratory glass bottle for dispersion of the immersion medium.

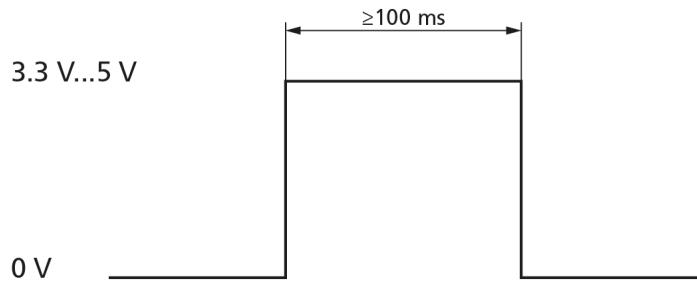Signal	Action
High	Switch the diaphragm pump on
Low	Switch the diaphragm pump off

Signal at the output **FEEDBACK_OUT**

The signal *FEEDBACK_OUT* provides information about whether immersion medium has been dispensed. An impulse is outputted for every time unit (0.1 s or 1.0 s) expired. This signal can be used to check whether immersion medium was dispensed in accordance with the required dispensing duration following transmission of one or more impulses to the input *DROP_IN*.

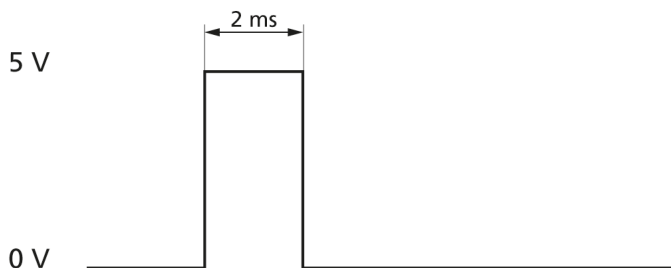Signal	Action
High	Immersion medium has been dispensed
Low	No immersion medium has been dispensed

7.5.3 Pin Assignment 5-Pin Port

Pin	Designation	Description	Function
1	DROP_IN	TTL-input (4.7 k Ω , 100 nF)	Dispense immersion medium
2	PUMP_IN	TTL-input (4.7 k Ω , 100 nF)	Switch diaphragm pump on/off
3	FEEDBACK_OUT	TTL-output (5 V)	Evaluate immersion medium dispensing
4	GND	GND	Ground
5	+5V	+5 V	Operating voltage
Housing	–	GND	Shielding

7.6 Operation at a PC (Instruction Set)

All the Liquid Dispenser settings can be configured and immersion medium dispensed by instruction set (ASCII instructions) via the USB interface.

You will find the instruction set on the installation CD:
» *Instruction Set Liquid Dispenser*

8 Cleaning

WARNING

Danger to life through the spread of substances which are harmful to health.

The Liquid Dispenser is not protected against corrosive, infectious, toxic, radioactive or other substances which are harmful to health.

- ▶ In the case of possible contamination with substances which are harmful to health, have the Liquid Dispenser professionally disinfected/decontaminated.

CAUTION

Danger of cutting by damaged laboratory glass bottles.

The laboratory bottles are made of glass and can be dropped and broken if not handled carefully. In this case, the immersion medium can leak.

- ▶ Make sure you handle the laboratory glass bottles carefully.

NOTICE

Possible equipment damage due to improper cleaning.

If the Liquid Dispenser is not sufficiently protected from moisture and liquids, this can lead to equipment damage.

- ▶ Disconnect the Liquid Dispenser from the power supply before starting cleaning work.
- ▶ Protect the ports and operating elements of the Liquid Dispenser from liquids.
- ▶ Ensure that no liquids enter the interior of the Liquid Dispenser.

8.1 Cleaning the Housing

If the Liquid Dispenser is dusty or soiled, clean the housing.

► Cleaning the housing

- 1 Disconnect the Liquid Dispenser from the power supply.
- 2 Clean the Liquid Dispenser with a lint-free cloth. With greater contamination, slightly dampen the cleaning cloth with water or a mild solvent-free disinfectant/cleaning agent.
- 3 Allow the Liquid Dispenser to dry.
- 4 Connect the Liquid Dispenser to the power supply.

Housing cleaning is now complete.

8.2 Removing the Immersion Medium

Remove the immersion medium if you wish to replace the immersion medium, clean all the parts containing immersion medium or not to use the Liquid Dispenser for a longer period of time, or to dispose of or dispatch it.

► Removing the immersion medium

- 1 Disconnect the Liquid Dispenser from the power supply.
- 2 Remove the bottle multiple distributors from the laboratory glass bottles 1 and 2 as follows: Hold the hose connections tight with one hand while turning the bottle multiple distributor counterclockwise with the other hand.
- 3 Dispose of the immersion medium from laboratory glass bottles 1 and 2 in accordance with manufacturer's instructions.
- 4 Fill laboratory glass bottle 1 up to the maximum filling height with a cleaning liquid suitable for the immersion medium.

- 5 Fix the bottle multiple distributors to laboratory glass bottles 1 and 2 as follows: Hold the hose connections tight with one hand while turning the bottle multiple distributor clockwise with the other hand.
- 6 Connect the Liquid Dispenser to the power supply.
- 7 Dispense cleaning liquid until the whole quantity of cleaning liquid is in laboratory glass bottle 2.
- 8 Dismantle the objective adapter from the objective.
» *Chapter 8.3 Dismantling the Objective Adapter, page 133*
- 9 Remove the bottle multiple distributors from the laboratory glass bottles 1 and 2.
- 10 Remove both laboratory glass bottles from the Liquid Dispenser.
- 11 Remove the cleaning liquid properly from the laboratory glass bottles.
- 12 Allow the laboratory glass bottles, the objective adapter and all hoses to dry.
- 13 Equip the Liquid Dispenser with the two laboratory glass bottles.
- 14 Fix the bottle multiple distributors to laboratory glass bottles 1 and 2.

The immersion medium has now been completely removed from the Liquid Dispenser.

The Liquid Dispenser can now be put back into service if required.


» *Chapter 5 Putting Into Service, page 92*

8.3 Dismantling the Objective Adapter

► Dismantling the objective adapter

- 1 Disconnect the Liquid Dispenser from the power supply.
- 2 Remove the immersion medium.
 - » *Chapter 8.2 Removing the Immersion Medium, page 131*
- 3 Open both bottle multiple distributors to allow the pressure to escape from the laboratory glass bottles.
- 4 Loosen the 4 screws ① to ④ evenly crosswise using the TORX wrench T8 without completely unscrewing the screws.

Notice: Never screw the screws out of the objective adapter completely since otherwise the seals may move under the screws.

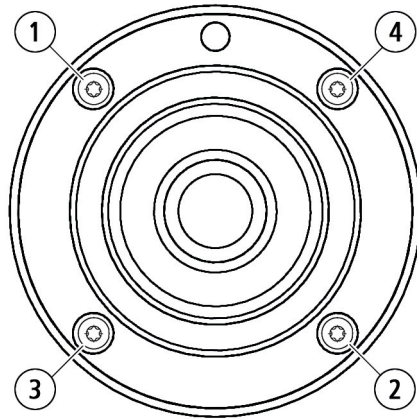- 5 Remove the objective adapter from the objective.
- 6 Remove the 2 hoses from the bottom of the objective adapter.
 - » *Chapter 3.3 Product Structure, page 82*

The objective adapter has now been completely dismantled.

9 Maintenance

Liquid Dispenser inverted for oil immersion

Replace/remove the immersion medium if

- the immersion medium is soiled.
- the specification of your immersion medium prescribes regular cleaning of the application system.
- the immersion medium you use leaves residue in the laboratory glass bottles or hoses.
- the Liquid Dispenser is not to be used for a longer period.
- the Liquid Dispenser is to be transported/disposed of.

Liquid Dispenser inverted for water immersion

Replace/remove the immersion medium if

- the immersion medium is soiled.
- the Liquid Dispenser is not to be used for a longer period.
- the Liquid Dispenser is to be transported/disposed of.

How to remove the immersion medium from the Liquid Dispenser is explained in the *Cleaning* chapter.

» *Chapter 8 Cleaning, page 130*

10 Troubleshooting

NOTICE

Possible equipment damage due to unauthorised repair attempts.

The Liquid Dispenser has not been designed for repair by the user. In case of violation, all warranty and liability claims will become void.

- ▶ Never open the housing of the Liquid Dispenser.
- ▶ Only have repair work carried out by our Customer Service Center.

The following overview can help you with remedying faults and malfunctions of the Liquid Dispenser.

If a fault or malfunction cannot be remedied using the overview, contact our Customer Service Center.

» *Chapter 11 Service, page 138*

Fault	Possible causes
Liquid Dispenser does not build up pressure	No power supply
	Hoses connected incorrectly
	Bottle multiple distributor fitted incorrectly
	Hoses damaged or leaking
	Liquid Dispenser faulty
No or continuous delivery of the immersion medium	Liquid Dispenser configured incorrectly
	Cabling faulty
	Hoses connected incorrectly
	Solenoid valve faulty
Start/stop button does not work, LED of the start/stop button does not light up/flash	Start/stop button configured incorrectly
	Start/stop button faulty

Remedy

Make sure the required operating voltage is connected.

» *Chapter 5.4 Supplying the Liquid Dispenser with Voltage, page 100*

Check the correct assignment of hoses and ports.

» *Chapter 3.3 Product Structure, page 82*

Check the correct assignment of bottle multiple distributors and laboratory glass bottles.

» *Chapter 3.3 Product Structure, page 82*

Contact our Customer Service Center.

» *Chapter 11 Service, page 138*

Contact our Customer Service Center.

» *Chapter 11 Service, page 138*

Configure the Liquid Dispenser again.

» *Chapter 6 Configuration, page 105*

- 1 Make sure that the data/signal cables used are shielded cables.
 - 2 If an external controller is being used, make sure the pin assignment is correct.
 - » *Chapter 7.5.3 Pin Assignment 5-Pin Port, page 129*
-

Check the correct assignment of hoses and ports.

» *Chapter 3.3 Product Structure, page 82*

Contact our Customer Service Center.

» *Chapter 11 Service, page 138*

Check the configuration of the start/stop button.

» *Chapter 6.1.3 Area 2 – Settings, page 108*

Contact our Customer Service Center.

» *Chapter 11 Service, page 138*

11 Service

The staff at our Customer Service Center will be happy to help with support queries, complaints and suggestions about our products and services.

Overview of the services from our Customer Service Center:

- Advice and support
- Support for customer-specific application questions
- Support for putting equipment into service
- Inspection and repair of equipment
- Fault analysis (also online via web conference)
- Delivery of spare parts and replacements
- User training (also online via web conference)

You can reach our Customer Service Center from 8 am to 4 pm, Monday to Friday.

Tel.: **+49 6441 9116-36**
 Fax: **+49 6441 9116-40**
 E-mail: **service@marzhauser.com**

So that we can help you quickly when you need service, please keep the following information about your Liquid Dispenser handy:

Description	Information on the nameplate
Model	
Art. No.	
S/N	

11.1 Replacement Components

Should you need to replace a component of the Liquid Dispenser at any time, please contact our Customer Service Center.

» *Chapter 11 Service, page 138*

Replacement components	Order no.
Laboratory glass bottle 100 ml	00-74-300-6810
Bottle multiple distributor 2×8	00-74-300-6812
Bottle multiple distributor 3×8	00-74-300-6813
USB cable	00-76-100-2802
Operating manual	00-74-300-9801
Installation CD	00-74-300-9802

12 Taking Out of Service


► Taking the Liquid Dispenser out of service

1 Clean the Liquid Dispenser completely.

» *Chapter 8 Cleaning, page 130*

Taking out of service is now complete.

13 Disposal

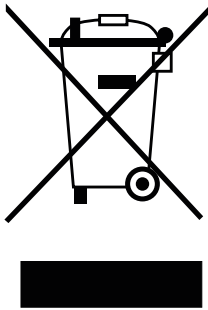The adjacent symbol means that the Liquid Dispenser is classified as electrical or electronic equipment in accordance with the EU Directive 2012/19/EU (WEEE) and thus cannot be disposed of with household waste.

Since the equipment is used commercially, it may not be disposed of via a communal collection point for electrical and electronic equipment.

Return the Liquid Dispenser to our Customer Service Center for proper disposal free of charge.

► Disposing of the Liquid Dispenser through Märzhäuser Wetzlar

- 1 Take the Liquid Dispenser out of service.
» *Chapter 12 Taking Out of Service, page 139*
- 2 Send the Liquid Dispenser for disposal to:
Märzhäuser Wetzlar GmbH & Co. KG
In der Murch 15
35579 Wetzlar (Germany)

All equipment purchased from Märzhäuser Wetzlar can be returned to our Customer Service Center at any time for free disposal. If you wish to dispose of the equipment yourself, always heed local laws and regulations.

14 Technical Data

Features

Operating element/ status display	Start/stop button with LED
Delivery principle	Pressure delivery
Operating pressure	< 500 mbar
Dosage	Time-controlled
Liquid quantity/flow	Programmable
Liquid return	Yes
Liquid reserve	100 ml
Operating voltage	5 V ($\pm 10\%$)
Power consumption	< 500 mA
Operation	Manual, via PC (software, ASCII instruction set), TANGO controller or external controller (TTL-signal transmission)
Interfaces	USB 2.0 type B (PC), 5-pin port (TANGO controller or external controller), 7-pin port (remote control)
Materials	Housing: Aluminium Hoses: PTFE, silicone Bottles: Glass
Protection type	IP41

Features

Dimensions (L × W × H)	Approx. 136 × 76 × 190 mm (approximate value incl. bottles and hoses, without supply lines to the objective adapter)
Weight	Approx. 1.3 kg (without immersion medium)

14.1 Power Adapter Recommendation

WARNING

Fire hazard and possible equipment damage caused by connection to a non-specified power supply.

The Liquid Dispenser has been designed for a power supply via USB or the 5-pin port. If the specifications for the power supply are not heeded, fire hazard and equipment damage can be the result.

- ▶ Make sure that your power supply complies with the specifications of the Liquid Dispenser.
 - » *Chapter 14 Technical Data, page 141*

For power supply you can use

- a PC with USB interface (USB 2.0 type A),
- a TANGO controller or an external controller with limited output current of max. 0.5...1 A or
- a 5 V USB power adapter with USB interface (USB 2.0 type A) and an output current of max. 0.5...1 A.

15 Index

A

- Accessories 78
- Area of application 79

C

- Cleaning 128
- COM port
 - Identifying* 115
 - Software settings* 105
- Customer Service Center 136

D

- Delivery contents 77
- Disposal 138

F

- Factory settings 117
- Faults 133
- Firmware 105

H

- Hose connections
 - Overview* 82
 - Setting up* 91

I

- Immersion medium
 - Filling* 96
 - Refilling* 118
 - Removing* 129
 - Sort* 84
- Installation requirements 85

L

- LED 81
 - configure* 107, 111
 - faulty* 134
 - flashes* 118

M

- Maintenance 132

N

- Nameplate 79

O

- Objective adapter
 - Dismantling* 131
 - Fitting* 94
 - Identifying* 79
- Operating elements 80
- Operating mode
 - Identifying* 118
 - Interval mode* 121
 - Manual mode* 119
- Operation 117

P

- Ports 80
- Power adapter recommend. 140
- Power supply 140
- Product structure 80
- Putting into service 90

R

- Repair 136
- Replacement components 137

S

- Safety 84
- Service 136
- Software
 - Config. by instruction set* 116
 - Installation* 102
 - Installing the driver* 100
- Start/stop button 81
 - configure* 107, 111
 - faulty* 134
- Switch-on 117

T

- Taking out of service 137
- Technical data 139
- Troubleshooting 133

U

- USB 82, 98, 140

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In der Murch 15
35579 Wetzlar (Germany)

Tel.: +49 6441 9116-0

Fax: +49 6441 9116-40

info@marzhauser.com

www.marzhauser.com

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© 2022 | All rights reserved. Any complete or partial copy requires written authorization from Märzhäuser Wetzlar GmbH & Co. KG.
All quoted product names are trademarks and/or registered trademarks of their respective owners.